

Operating Instructions

LAUDA Compact Thermostats
C 6 CP, C 12 CP, C20 CP
K 6 KP, K 12 KP, K 20 KP
LAUDA Clear-View Thermostats
D 15 KP, D 20 KP, D 30 KP,
DL 15 KP, DL 20 KP

| | | |
|----------|--|-----------|
| 1 | BRIEF OPERATING INTRUCTIONS | 6 |
| 1.1 | Fitting the tubing to the pump connections:..... | 6 |
| 1.2 | Operating safety | 7 |
| 1.3 | Warning | 7 |
| 2 | TECHNICAL DATA (TO DIN 58966) | 8 |
| 3 | GENERAL CONSTRUCTION AND TECHNICAL DESCRIPTION..... | 13 |
| 3.1 | Operating principle | 13 |
| 3.1.1 | Bath/Circulation Thermostats..... | 13 |
| 3.1.2 | Clear-View Thermostats | 13 |
| 3.2 | Materials..... | 13 |
| 3.3 | Cooling coil | 14 |
| 3.4 | Pumps | 14 |
| 3.5 | Temperature control and electronics | 15 |
| 3.6 | Mains supply output 34 H..... | 15 |
| 3.7 | Controlled cooling..... | 15 |
| 3.8 | Remote operation (FBC) (option) | 15 |
| 4 | SAFETY DEVICES AND WARNING NOTES | 16 |
| 4.1 | Safety functions | 16 |
| 4.2 | Why can a thermostat be dangerous? | 16 |
| 4.3 | Important notes | 17 |
| 4.4 | Warning notes | 18 |
| 4.4.1 | Temperatures..... | 18 |
| 4.4.2 | Mains connection | 18 |
| 4.4.3 | Mains cable | 18 |
| 4.4.4 | Fume extraction | 18 |
| 4.4.5 | Cooling water, steam production..... | 18 |
| 5 | BATH LIQUIDS AND HOSE CONNECTIONS | 19 |
| 5.1 | Bath liquids..... | 19 |
| 5.2 | Hose connections | 20 |
| 6 | UNPACKING, ASSEMBLY AND SETTING UP..... | 21 |
| 6.1 | Unpacking | 21 |
| 6.2 | Setting up, operation as bath thermostat | 21 |
| 7 | CONNECTION OF EXTERNAL SYSTEMS..... | 22 |
| 7.1 | Closed external circuits..... | 22 |
| 7.2 | Open systems (baths)..... | 22 |
| 7.2.1 | Fitting the level controller | 23 |
| 8 | COOLING THE THERMOSTATS..... | 24 |

| | | |
|----------|--|-----------|
| 8.1 | Mains water cooling..... | 24 |
| 8.2 | Through-flow chillers DLK 10, DLK 25 and DLK 45..... | 24 |
| 9 | STARTING UP | 25 |
| 9.1 | Filling..... | 25 |
| 9.2 | Connection to supply..... | 25 |
| 9.3 | Basic functions..... | 25 |
| 9.3.1 | Supply switch-on | 25 |
| 9.3.2 | Standard display..... | 26 |
| 9.3.3 | Basic action on inputs and outputs | 26 |
| 9.3.4 | Overtemperature switch-off point..... | 27 |
| 9.3.5 | Low temperature switch-off point | 28 |
| 9.3.6 | Setpoint input | 28 |
| 9.4 | Parameter level PAR | 29 |
| 9.4.1 | Auto-adaptation..... | 29 |
| 9.4.2 | Output limitation | 29 |
| 9.4.3 | Display resolution L1 | 30 |
| 9.4.4 | Contact input Fault 14 N | 30 |
| 9.4.5 | Baudrate RS 232..... | 31 |
| 9.4.6 | Menu language | 31 |
| 9.4.7 | Calibrating the analogue output channels..... | 32 |
| 9.4.8 | Operation with through-flow chiller DLK 45 with proportional cooling | 32 |
| 9.5 | Calibration of the temperature measurement circuits..... | 33 |
| 9.6 | Control parameters | 34 |
| 9.6.1 | Indication and input of the control parameters..... | 34 |
| 9.6.2 | Recommendations for the control parameters..... | 34 |
| 9.6.3 | Bath temperature limitation | 35 |
| 9.6.4 | Correction limitation | 35 |
| 9.7 | External control | 36 |
| 9.7.1 | External measurement inputs and external controller..... | 36 |
| 9.7.2 | Start of external control | 37 |
| 9.7.3 | Notes..... | 37 |
| 9.8 | Working with controlled cooling | 38 |
| 9.9 | Operation with programmer | 40 |
| 9.9.1 | Programme input | 40 |
| 9.9.2 | Example of a programme | 42 |
| 9.9.3 | Programm test | 43 |
| 9.9.4 | Changing the programme data | 43 |
| 9.9.5 | Programme start, interruption and abort..... | 43 |
| 9.10 | Connection for analogue signals socket 52 S..... | 43 |
| 9.11 | Analogue inputs | 44 |
| 9.12 | Analogue outputs..... | 48 |
| 9.12.1 | Temperature signal channel 1 | 48 |
| 9.12.2 | Temperature signal channel 2 | 50 |
| 9.13 | Safety function | 50 |
| 9.13.1 | Low-level cut-out | 50 |
| 9.13.2 | Adjustable overtemperature limiter | 50 |
| 9.13.3 | Connection potential-free contact "Combined fault" 12 N (Alarm off)..... | 52 |
| 9.14 | Serial interface RS 232 C | 52 |
| 9.14.1 | Data of the RS 232 C interface | 52 |

| | | |
|-------------|---|-----------|
| 9.14.2 | General principles | 53 |
| 9.14.3 | Output commands..... | 53 |
| 9.14.4 | Requesting data from the thermostat..... | 55 |
| 9.14.5 | Error messages on the computer | 57 |
| 9.15 | Timing clock function | 57 |
| 9.15.1 | Setting and indication of date and time | 57 |
| 9.15.2 | Timing clock function | 58 |
| 9.15.3 | Activating the timing clock function | 58 |
| 10 | MAINTENANCE | 60 |
| 10.1 | Safety notes in case of repairs | 60 |
| 10.2 | Repair and re-initialisation | 60 |
| 10.3 | Cleaning | 61 |
| 10.4 | Spares ordering..... | 61 |
| 11 | ACCESSORIES | 62 |

1 Brief operating instructions

Even if you find these brief instructions initially sufficient please read the following sections, especially Section 4: "Safety devices and warning notes".

For safe operation of the equipment it is essential that the information in these Operating Instructions is observed.

Check the thermostat and the accessories during unpacking for any transport damage and if necessary inform the carrier or the postal authority.

Assemble the unit according to Section 6 and add extra items as appropriate.

1.1 *Fitting the tubing to the pump connections:*

Without external system: for improved circulation within the bath remove the closing plugs from both of the pump connections, fit the tubing nipples and link them together with e.g. Perbunan tubing (up to 120°C) or better a metal tubing.

For clear-view thermostats: for improved circulation within the bath close the pump connections with the closing plugs (see "Standard accessories", Item 6.1).

With external system: make tubing connections to the external system.

Protect the tubing with hose clips against slipping off.



When working near the ambient temperature connect up the external cooling according to Section 8.

Use only softened water or LAUDA bath liquids (Section 5). Fill up the bath to a level about 2 cm below the cover plate.

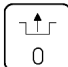
Check the supply voltage against the details on the label. Insert the mains plug.

Switch on the supply switch (green lamp lights up)

The display shows the software version and the type of unit, followed by the standard display.


Select the required indications using the keys  and  in the SHIFT mode. It is useful to show the setpoint (Ts) in display line 2 (L2) (see Items 9.3.1 - 9.3.3).

Set the overtemperature switch-off point (TO) slightly above the operating temperature.

If there is an error message, press the key  and perhaps increase TO.

Important: TO must be at least 25 K below the fire point of the bath liquid used (see Item 9.3.4)!

When connecting up an external system, check that filling this system does not cause the level inside the thermostat to fall more than is permitted.

When the thermostating liquid has reached the setpoint the symbol  starts to flash on the right in line 1 (L1) of the display.

After the unit has stabilised the bath temperature (T_i) corresponds to the setpoint (T_s).

1.2 *Operating safety*

The thermostat must be operated only with non-flammable bath liquids or with flammable bath liquids up to 25°C below their fire point, otherwise there is the possibility that a flammable atmosphere may form (see Item 4.2).

1.3 *Warning*

Parts of the bath cover may reach temperatures above 70°C when working at higher temperatures. The outflow and return pipes of the pumps reach the operating temperature. Touching them is dangerous because of very high or low temperatures!

2 Technical Data (to DIN 58966)

| Type | | C 6 CP | C 12 CP | C 20 CP | K 6 KP | K 12 KP | K 20 KP |
|--|------------------|--|--------------------|--------------------|-------------------|----------------|-------------------|
| Operating temperature range | (°C) | -30...300 | | | | | |
| Working temperature range with water cooling (water 15°C) | (°C) | 45...300 | 40...300 | 35...300 | 50...300 | 40...300 | 45...300 |
| | (°C) | 20...300 (> item 4.4.5) | | | | | |
| Ambient temperature range | (°C) | 5...40 | 5...40 | 5...40 | 5...40 | 5...40 | 5...40 |
| Temperature setting / Resolution | | membrane keypad with 16 keys, setpoint input with 0.01°C resolution | | | | | |
| Bath temperature measurement | | built-in digital thermometer with 0.01°C resolution; accuracy and stability of measurement (electronics without probe) better than 0.05% ± 0.05 K*). Temperature probe Pt 100 to DIN IEC 751, Class B, can be calibrated additively at each measurement point | | | | | |
| Display | | back-lit LCD matrix display, 2 lines with 16 signs each, 10 mm character height | | | | | |
| External temperature measurement | | 2 separate temperature measurement circuits for external Pt 100 to DIN IEC 751 in 4-wire circuit, accuracy and stability of measurement better than 0.05% ± 0.05 K*), can be calibrated at each measurement point. One probe produces the measurement for external control | | | | | |
| Temperature control | | modified PID controller with automatic structure selection. Control parameters through auto-adaptation or manual input with external control a cascade controller operates on the measurement of one of the two external circuits (T1 or T2). | | | | | |
| Temperature variation at 70°C in the bath | (°C) | 0.01*) | | | | | |
| Heater power | (kW) | 2 | 2 | 2 | 3 | 3 | 3 |
| Effective surface area of cooling coil | (cm²) | 350 | | | | | |
| Safety system | | Overtemperature protection with selctable switch off oint, low-level protection to EN 61010 FL (DIN 12879 Class 2) | | | | | |
| Pump output against zero head Pressure/Suction max. pump pressure | (L/min) (bar) | 20/16 0.32/0.25 | 20/16 0.32/0.25 | 20/16 0.32/0.25 | 24/18 0.5/0.34 | 24/- 0.5/- | 24/18 0.5/0.34 |
| Pump connections | (I.W.) | M16x1; nipples 13 Ø | | | | | |
| Filling volume max. | (L) | 4...6 | 9...14 | 14...21 | 4.5...7.5 | 6...12 | 11...18 |
| Bath opening (B x T) | (mm) | 150x130 | 300x175 | 300x350 | 150x130 | 150x130 | 300x175 |
| Bath depth | (mm) | 160 | 160 | 160 | 200 | 320 | 200 |
| Usable liquid depth | (mm) | 140 | 140 | 140 | 180 | 300 | 180 |
| Height to top of bath | (mm) | 220 | 220 | 220 | 260 | 380 | 260 |
| Overall size (W x D x H) | (mm) | 200x350x420 | 375x415x420 | 375x590x420 | 200x350x460 | 225x375x580 | 375x415x460 |
| Weight | (Kg) | 12 | 19 | 23 | 14 | 19 | 22 |
| Supply | (V;Hz) | 230;50 / 230;60 protection class 1 to VDE 0106 | | | | | |
| Loading (max.) without accessories | (kW) | 2.2 | 2.2 | 2.2 | 3.2 | 3.2 | 3.2 |
| Current take up (max.) with accessories | (A) | 15 | 15 | 15 | 16 | 16 | 16 |
| Nominal current of fuses on customer's side | (A) | T 16 A | | | | | |
| | | Units are conform to EU Guideline 89/336/EWG (EMC) and 73/23/EWG (low-voltage) and carry the CE mark (230 V; 50 Hz) | | | | | |
| Cat-No.: | | | | | | | |
| 230V; 50Hz | | LCB 156 | LCB 158 | LCB 160 | LCB 164 | LCB 166 | LCB 168 |
| 230V; 60Hz | | LCB 256 | LCB 258 | LCB 260 | LCB 264 | LCB 266 | LCB 268 |

Units of different power supplies have different heating capacities as well as different values for power consumption (see type label)!

Technical changes reserved!

*) > item 4.3

| Type | | D 15 KP | D 20 KP | D 30 KP | DL 15 KP | DL 20 KP |
|--|------------------|--|----------------------------|----------------------------|----------------------------|----------------------------|
| Operating temperature range | (°C) | 0...230 | 0...230 | 0...230 | -40...100 | -40...100 |
| Working temperature range | (°C) | 45...230 | 40...230 | 40...230 | 45...100 | 45...100 |
| with water cooling (water 15°C) | (°C) | 20...300 (≥ item 4.4.5) | 20...300 (≥ item 4.4.5) | 20...300 (≥ item 4.4.5) | 20...100 (≥ item 4.4.5) | 20...100 (≥ item 4.4.5) |
| Ambient temperature range | (°C) | 5...40 | 5...40 | 5...40 | 5...40 | 5...40 |
| Temperature setting / Resolution | | membrane keypad with 16 keys, setpoint input with 0.01°C resolution | | | | |
| Bath temperature measurement | | built-in digital thermometer with 0.01°C resolution; accuracy and stability of measurement (electronics without probe) better than 0.05% ± 0.05 K*). Temperature probe Pt 100 to DIN IEC 751, Class B, can be calibrated additively at each measurement point | | | | |
| Display | | back-lit LCD matrix display, 2 lines with 16 signs each, 10 mm character height | | | | |
| External temperature measurement | | 2 separate temperature measurement circuits for external Pt 100 to DIN IEC 751 in 4-wire circuit, accuracy and stability of measurement better than 0.05% ± 0.05 K*), can be calibrated at each measurement point. One probe produces the measurement for external control | | | | |
| Temperature control | | modified PID controller with automatic structure selection. Control parameters through auto-adaptation or manual input with external control a cascade controller operates on the measurement of one of the two external circuits (T1 or T2). | | | | |
| Temperature variation at 70°C in the bath | (°C) | 0.01*) | | | | |
| Heater power | (kW) | 3 | 3 | 3 | 2 | 2 |
| Effective surface area of cooling coil | (cm²) | 350 | | | | |
| Safety system | | Overtemperature protection with selctable switch off oint, low-level protection to EN 61010 FL (DIN 12879 Class 2) | | | | |
| Pump output against zero head Pressure/Suction max. pump pressure | (L/min) (bar) | 18 0.4 | 18 0.4 | 15 0.3 | 18 0.4 | 18 0.4 |
| Pump connections | (I.W.) | M 16x1; Oliven 13Ø | | | | |
| Filling volume max. | (L) | 13...15 | 21...24 | 32...36 | 13...15 | 21...24 |
| Bath opening (B x T) | (mm) | 230x135 | 405x135 | 585x135 | 230x135 | 405x135 |
| Bath depth | (mm) | 320 | 320 | 320 | 320 | 320 |
| Usable liquid depth | (mm) | 265 | 265 | 265 | 265 | 265 |
| Height to top of bath | (mm) | 390 | 390 | 390 | 395 | 395 |
| Overall size (W x D x H) | (mm) | 480x225x590 | 710x235x590 | 1010x235x590 | 506x250x595 | 740x250x595 |
| Weight | (Kg) | 12 | 19 | 23 | 14 | 19 |
| Supply | (V;Hz) | 230;50 / 230;60 protection class 1 to VDE 0106 | | | | |
| Loading (max.) without accessories | (kW) | 3.2 | 3.2 | 3.2 | 2.2 | 2.2 |
| Current take up (max.) with accessories | (A) | 16 | 16 | 16 | 16 | 16 |
| Nominal current of fuses (max.) on customer's side | (A) | T 16 A | | | | |
| | | Units are conform to EU Guideline 89/336/EWG (EMC) and 73/23/EWG (low-voltage) and carry the CE mark (230 V; 50 Hz) | | | | |
| Cat-No: | | | | | | |
| 230V; 50Hz | | LCD 121 | LCD 122 | LCD 123 | LCD 1272 | LCD 1273 |
| 230V; 60Hz | | LCD 221 | LCD 222 | LCD 223 | LCD 2272 | LCD 2273 |

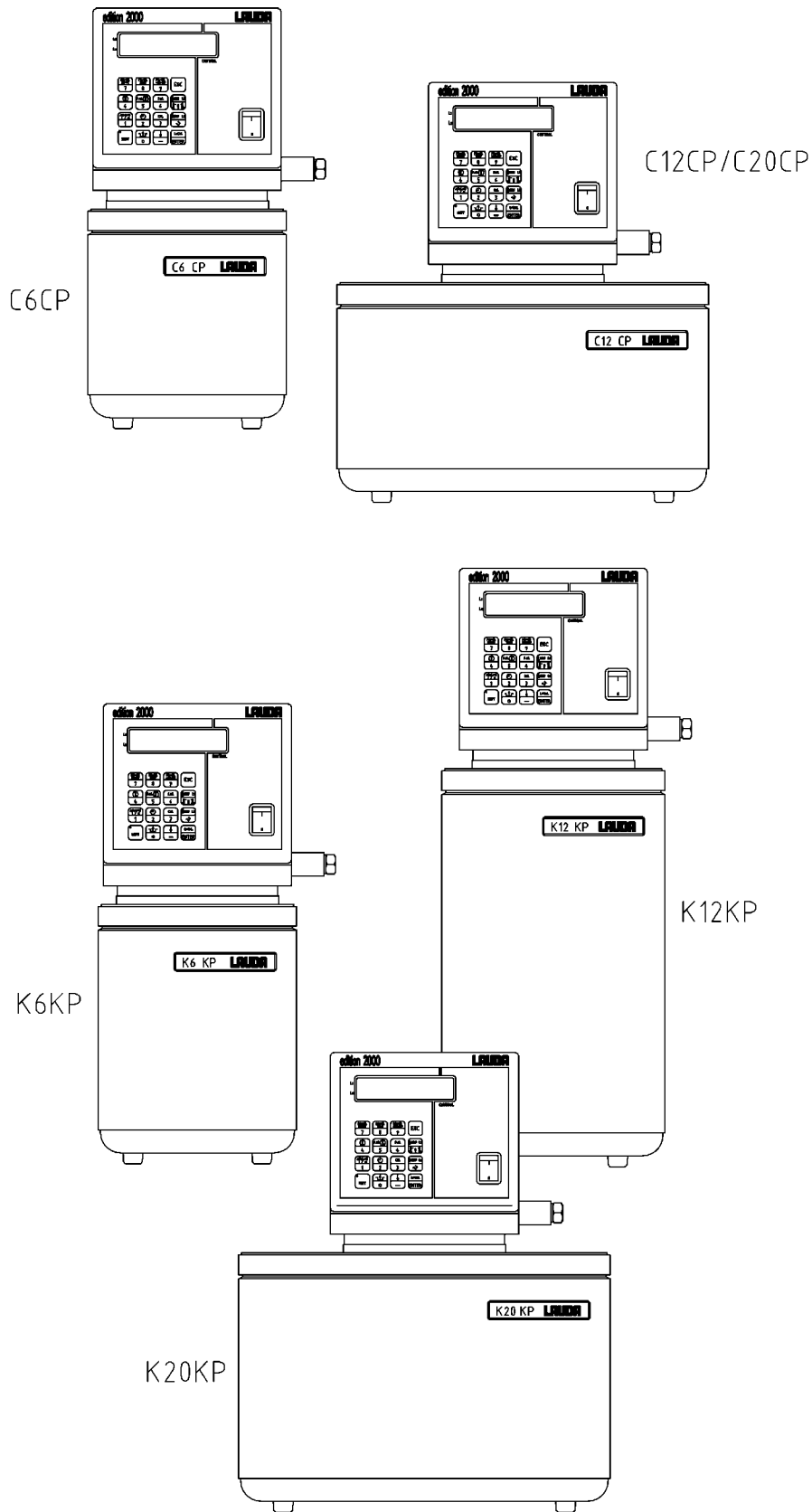
Units of different power supplies have different heating capacities as well as different values for power consumption (see type label)!

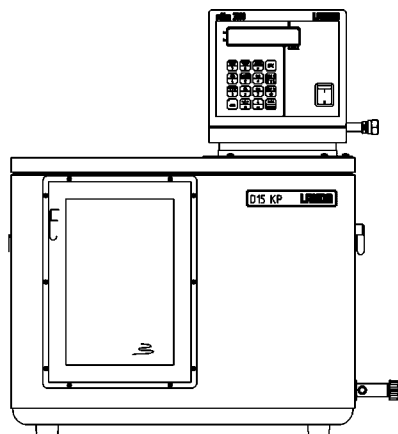
Technical changes reserved!

*) ≥ item 4.3

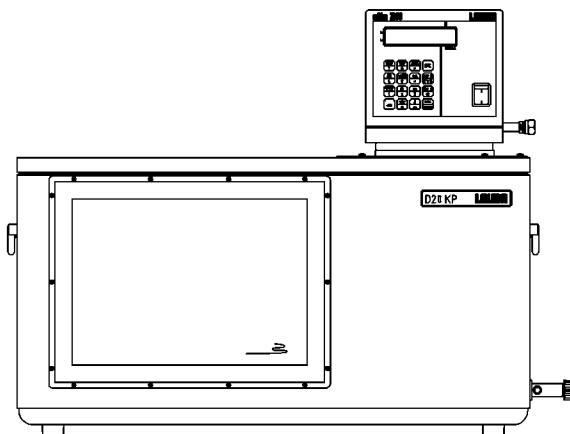
LAUDA Compact Thermostats

C 6 CP, C 12 CP, C 20 CP, K 6 KP, K 12 KP, K 20 KP,
D 15 KP, D 20 KP, D 30 KP, DL 15 KP, DL 20 KP

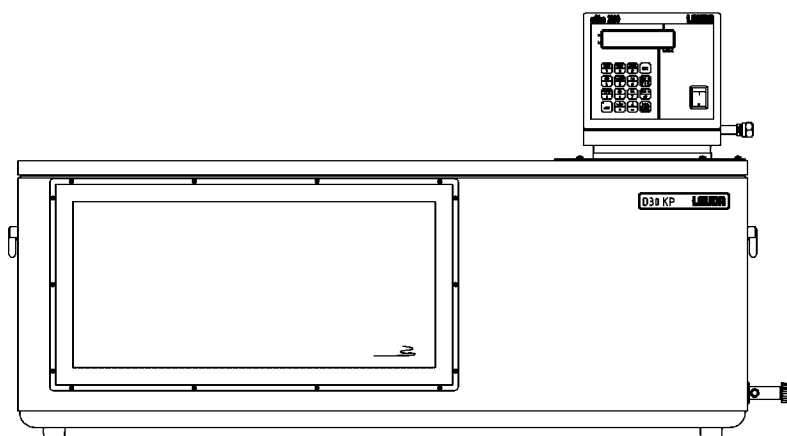




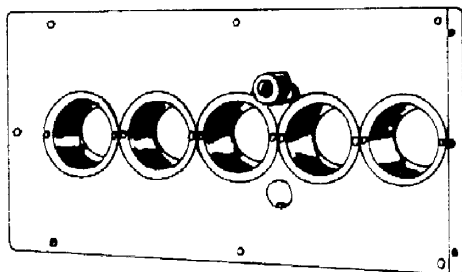
D15KP, DL15KP



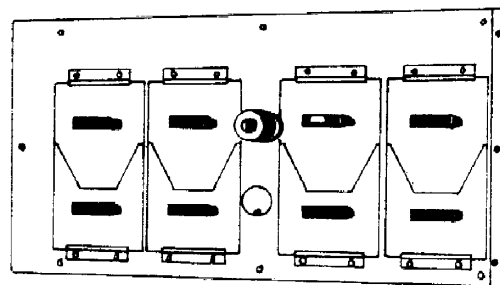
D20KP, DL20KP



D30KP

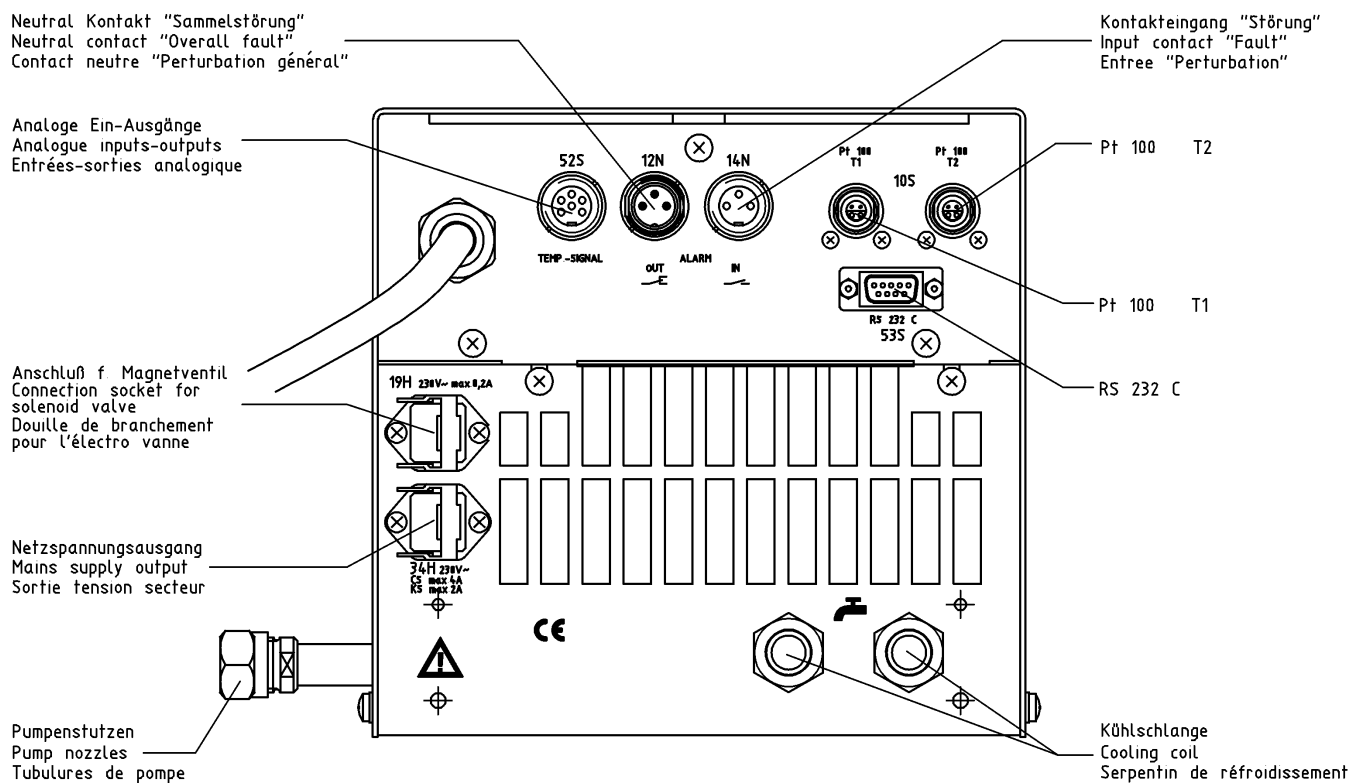
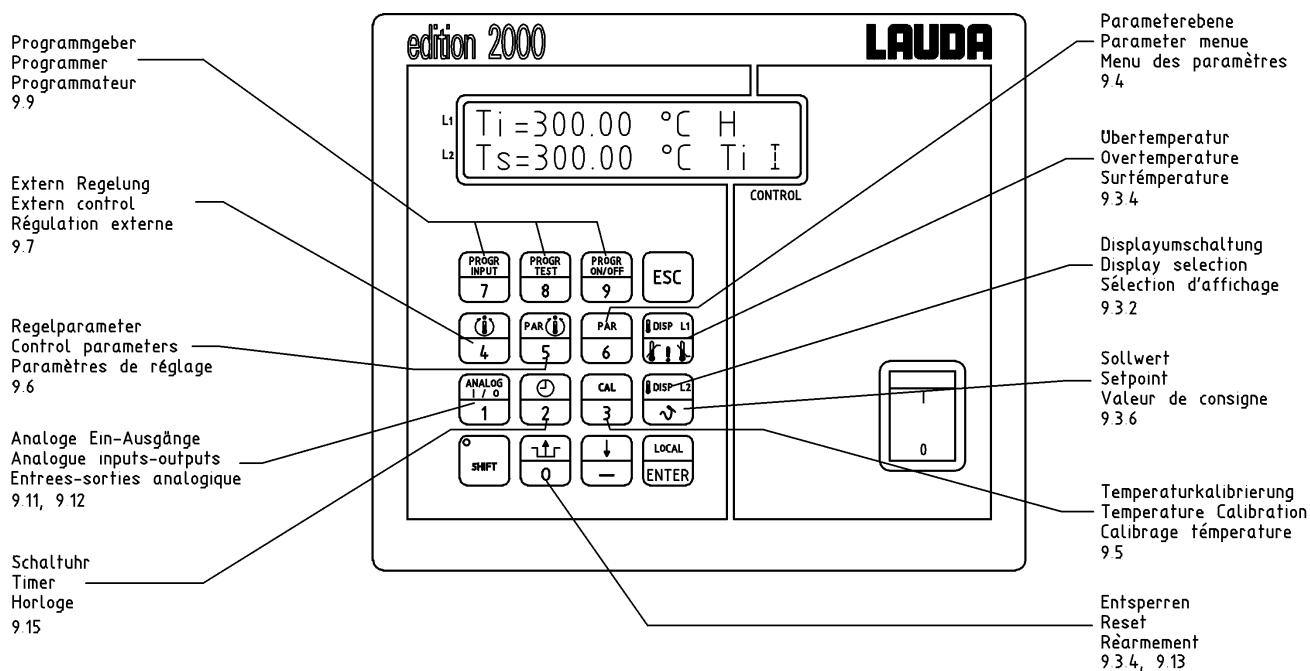


Deckplatte für Kapillarviskosimeter
 Cover plate for capillary viscosimeter
 Plaque support pour viscosimetres capillaires



Deckplatte für VISCOTIMER-Messtativ
 Cover plate for VISCOTIMER measuring stands
 Plaque support pour statifs de mesure du VISCOTIMER

LAUDA Compact Thermostats
C 6 CP, C 12 CP, C 20 CP, K 6 KP, K 12 KP, K 20 KP,
D 15 KP, D 20 KP, D 30 KP, DL 15 KP, DL 20 KP



3 General construction and technical description

3.1 Operating principle

3.1.1 Bath/Circulation Thermostats

The LAUDA bath/circulation thermostats Series C and K with Electronics P differ in bath volume, bath depth, pump type and output as well as in heating capacity. All types offer an operating temperature range of $-30\ldots300^{\circ}\text{C}$. Laboratory thermostats operate with liquids (operating medium, heat transfer fluid) which serve for energy transfer to the product to be thermostated.

The thermostated products can be immersed in the thermostatic bath (bath thermostat), or placed in an external open bath whose liquid is circulated by the pump of the thermostat.

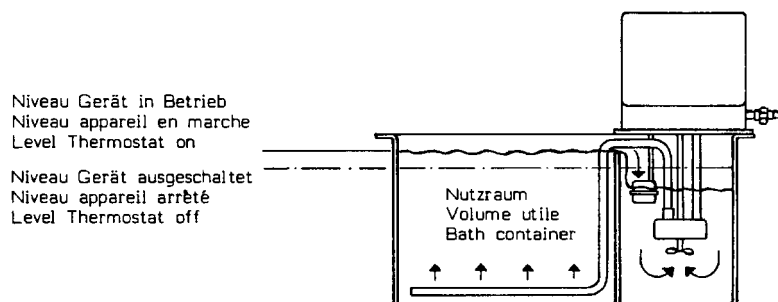
When operating as circulator the thermostatic liquid is pumped through an external heat exchanger arranged by the user in which a product is being thermostated (jacketed vessels, reactors, heat exchangers).

3.1.2 Clear-View Thermostats

The LAUDA clear-view thermostats Series D and DL with Electronics P differ in bath volume and temperature range. The units D...KP are equipped with double-walled inspection glasses and offer an operating temperature range of $0\ldots230^{\circ}\text{C}$. The units DL...KP are equipped with four-fold insulating glasses for working in an operating temperature range of $-40\ldots100^{\circ}\text{C}$.

The units are especially suitable for thermostating capillary viscometers.

The bath is divided into a bath container which is supplied with inspection glasses of several layers at the front and the rear, and an ante-chamber, in which the function elements (heater, pump, cooling coil etc.) are located. This provides a constant bath level in the bath container.



3.2 Materials

All materials in contact with the bath liquid are made from rust-free stainless steel or materials of similar anti-corrosion properties.

3.3 Cooling coil

All units are fitted with a cooling coil which permits cooling, e.g. with water, for working temperatures in the range of the ambient temperature (see Item 4.4.5 and Section 8).

3.4 Pumps

All units except Type K 12 KP and the clear-view thermostats are equipped with a centrifugal pressure/suction pump. This can be used to operate both external open baths as well as closed external systems (reactors).

The immersion pumps are supplied in the two performance classes C and K. The thermostat K 12 KP is fitted with a pure pressure pump since this unit operates specially in the high temperature range with closed external systems.

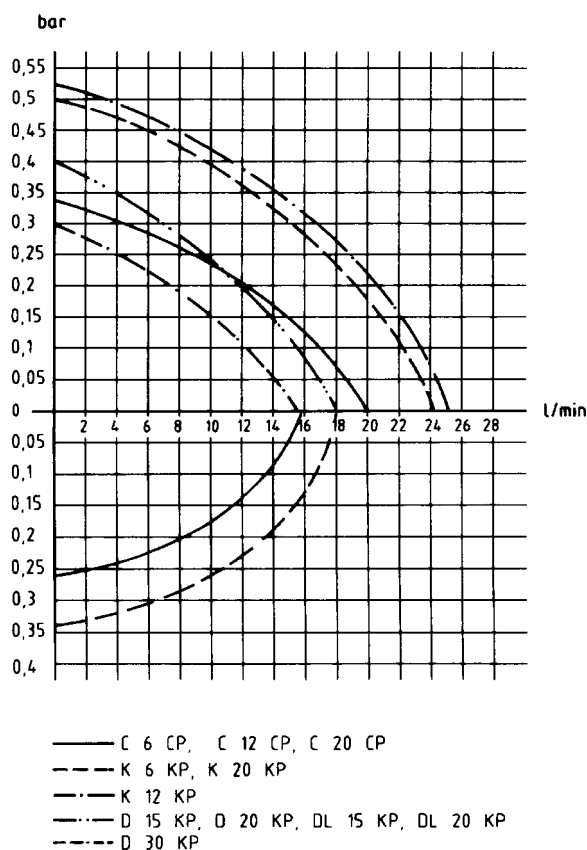
The pressure pumps of the clear-view thermostats ensure an excellent circulation using jet pipes and can also be used for external circulation if necessary.

The pumps are driven by external-runner motors with a continuous shaft.

The pumps operate perfectly up to a viscosity of approx. 70 mm²/sec. (K 12 KP approx. 120 mm²/sec.), with the pump output decreasing rapidly with increasing viscosity.

Pump characteristics

230V; 50Hz



3.5 *Temperature control and electronics*

The units operate with a Pt 100 resistance thermometer for measuring the bath temperature (Ti). The bath temperature, all other temperature values and message signals as well as inputs are indicated as 2 x 16 characters (10 mm high) on a liquid crystal display (LCD) with background illumination. Input of the setpoint (Ts) and of all other parameters is made by using a membrane keypad with 16 keys and the operator guidance in the LCD display field. All inputs are stored even when the thermostat is switched off or if the supply fails.

The digitizing of the Pt 100 resistance signal is performed in the microprocessor by continuous comparison with precision resistors. The secondary control using a modified PID control algorithm is purely digital. The tubular heater for the heating of the bath is then operated electronically using a triac with burst firing action. The tubular heaters have a surface loading of approx. 6W/cm².

3.6 *Mains supply output 34 H*

The 230V supply voltage is available at the socket 34 H at the back in normal operation and with the unit switched on. The maximum current which can be drawn there is 4A on C-P units and 2A on K-P units. In case of a fault this supply is switched off. This output can be used e.g. to connect a through-flow chiller or a non-return fitting (Cat. No. UD 125).

Suitable mating plug

Cat. No. EQS 045

3.7 *Controlled cooling*

The units are equipped for controlled cooling to operate a solenoid valve which controls the cooling water flow. This provides fully automatic cooling (20...100°C). It ensures faster heating up (compared to continuous cooling), greatly reduced water consumption and improved temperature control during heat dissipation since the heater does not operate against the cooling action.

Solenoid valve for cooling water control

Cat. No. UD 085

3.8 *Remote operation (FBC) (option)*

As an option the units can be converted to remote control; then the entire electronics with control panel is removed from the unit and used for remote operation. An adapter for the cable connections is required on the basic unit and the control panel is placed in an extra housing. The conversion is to be carried out by a qualified electrician only. All necessary components except the connection cables are supplied as part of the kit. Please specify the length of the connection cables. See accessories in the appendix.

4 Safety devices and warning notes

4.1 Safety functions

The built-in overtemperature limiter can be set over the complete operating temperature range.

The bath temperature is measured by a separate Pt 100 resistance sensor (T_{Si}) and processed by a separate analogue/digital converter. This measured value is compared with the measured value of the bath temperature probe (T_i) continuously. If the measurements differ by more than ± 15 K the thermostat switches off as in the case of a low-level or overtemperature fault.

The function of the microprocessor is monitored by an integrated watchdog circuit and an additional counter which operates similarly to a normal watchdog circuit but is also capable of switching off the unit in case of a strobe failure.

When the set overtemperature switch-off point (T_O) is exceeded the unit switches off permanently on all poles.

A float switch with magnetic coupling acts as a low-level cut-out and also switches off the unit (pump and heater) permanently on all poles.

In both fault conditions the display shows the corresponding message, and additionally an audible signal draws attention to the fault. The switch-off function of the safety circuit remains stored even during a break in the supply or after having switched off the supply.

Reset is possible by pressing the reset key , but only after having eliminated the troubles.

The pump motor is fitted with a temperature monitor which switches off if the motor winding overheats. The heater is also switched off simultaneously. After the motor winding has cooled down the pump starts up automatically.

4.2 Why can a thermostat be dangerous?

1. Thermostats are equipped with heaters supplying the necessary heat to the thermostating liquid. If the temperature control fails or if the liquid level is too low, the heater may reach temperatures which can lead to a fire in the laboratory, especially in combination with flammable liquids.
2. When using the thermostat as a circulation thermostat a hose may break, causing hot liquid to spill and endangering people and goods.

The safety requirements on thermostats therefore depend on whether

- o non-flammable or flammable liquids are used
- o operation is with or without supervision.

The thermostats described in these Operating Instructions are protected against overtemperature and low liquid level when operated according to the regulations (FL).

The units can be operated with non-flammable bath liquids and with flammable bath liquids up to 25°C below their fire point (EN 61010), on condition that there is a correct adjustment and regular testing (see Item 9.13) of overtemperature and low-level protection.

4.3 Important notes

The user is only protected against those hazards which are caused by exceeding the temperature and by low liquid level.

Further hazards may arise from the type of product being thermostated, e.g. a shift above or below certain temperature levels or a fracture of the container and a reaction with the thermostatic liquid etc.

It is impossible to cover all possible causes, and they remain largely within the decision and responsibility of the user.

Values for temperature variation and indication accuracy apply under normal conditions according to DIN 58966. In special cases high-frequency electromagnetic fields may lead to less favourable values. There is no loss of safety.

Units are in accordance with EMC directive EN 61 326-1, class A *:

Units are only suitable for use in industrial areas as disturbing voltage fluctuations might occur.

Note: The units must only be used according to the descriptions indicated in these Operating Instructions.

This includes operation by properly qualified and instructed personnel.

The units are not designed for operation under medical conditions according to EN 60601-1 or IEC 601-1!

4.4 *Warning notes*

4.4.1 Temperatures

Parts of the bath cover may reach temperatures above 70°C when working at higher temperatures. The outflow and return pipes of the pumps reach the operating temperature. Touching them is dangerous because of very high or low temperatures!

4.4.2 Mains connection

Connect the unit only to mains sockets with protective earth contact (PE) which must not have a fuse higher than T 16 A!

4.4.3 Mains cable

We have ensured that the mains cable and other plug connections do not touch any hot parts. Please check that there is no contact between the connecting tubings filled with hot liquid or other hot parts and the mains cable!

4.4.4 Fume extraction

Depending on the bath liquid used and the operating method there is a possibility that toxic vapours may be produced. In that case it is necessary to provide an appropriate fume extraction. Pull out the mains plug before cleaning the bath with solvents. Provide appropriate fume extraction. Before starting up the unit it is absolutely essential to ensure that the bath contains no explosive mixture. If necessary purge it with nitrogen!

4.4.5 Cooling water, steam production

Use cooling coils with cooling water only at operating temperatures below 100°C; at higher temperatures there is a danger that superheated steam may be produced. When changing the bath liquid from water to heat transfer fluids for temperatures above 100°C any remaining water - including the one in the hoses and external system - has to be removed completely. Otherwise there is a danger of burns because of delayed boiling.

* Notice only valid for EU countries

5 Bath liquids and hose connections

The operating temperature ranges of the bath liquids and hoses are for general information only and may be restricted through the operating temperature range of the units or the safety requirements specified in the appropriate standards (see Item 4.2).

5.1 Bath liquids

| LAUDA Designation | | Working-temperature-range | Chem. Designation | Viscosity (kin) | Viscosity (kin) at temperature | Fire-point | Ref.No. Quantity | | |
|-------------------|---------------------|---------------------------|----------------------------|--------------------|--------------------------------|------------|------------------|---------|---------|
| | Former designation | from °C to °C | at 20°C | mm ² /s | mm ² /s | | 5 l | 10 l | 20 l |
| | Water | +5...+90 | deionised water ① | -- | -- | -- | | | |
| Kryo 30 ② | G 100 ② | -30...+90 | Mono-ethylene-glycol/water | 4 | 50 at -25°C | -- | LZB 109 | LZB 209 | LZB 309 |
| Kryo 85 | Ultra-Therm XLT | -85...+30 | Silicone oil | 1.76 | 17 at -80°C | > 56 | LZB 113 | LZB 213 | LZB 313 |
| Ultra 350 ③ | 330 SCB | +30...+200 | synthetic thermal oil | 47 | 28 at +30°C | > 240 | LZB 107 | LZB 207 | LZB 307 |
| Ultra 300 | Ultra-Therm SW 300N | +80...+300 | Silicone oil | 170 | 35 at +80°C | > 400 | LZB 108 | LZB 208 | LZB 308 |
| Therm 180 | RDS 20 | 0...+180 | Silicone oil | 23 | 25 at 20°C | > 288 | LZB 114 | LZB 214 | LZB 314 |
| Therm 230 | RDS 50 | +60...+230 | Silicone oil | 54 | 28 at +60°C | > 362 | LZB 117 | LZB 217 | LZB 317 |




- ① At higher temperatures → Evaporation losses → Use bath covers (➤ Section 10. Accessories). Distilled water or fully deionised water must only be used with the addition of 0.1g sodium carbonate (Na₂CO₃) /l water, otherwise → danger of corrosion!
- ② Water content falls after prolonged operation at higher temperatures → mixture becomes flammable (flash point 128 °C). → Check the mixture ratio with a densitometer.
- ③ Do not use in conjunction with EPDM tubing!
 - When selecting bath liquids it should be noted that performance must be expected to worsen at the lower limit of the operating temperature range due to increasing viscosity. The full operating range should only be utilised if really necessary.
 - The operating ranges of the bath liquids and tubing represent general data which may be limited by the operating temperature range of the unit.



Silicone oil causes pronounced swelling of Silicone rubber → never use Silicone oil with Silicone tubing!

DIN Safety data sheets are available on request!

5.2 Hose connections

| Tubing type | Int. dia. Ø mm | Temperature range °C | Application | Cat. No. |
|---|--------------------------------|-------------------------|---|----------------|
| EPDM-tubing uninsulated | 9 | 10...120 | for all bath liquids except Ultra 350 and mineral oils | RKJ 111 |
| EPDM-tubing insulated | 9 | -60...120 | for all bath liquids except Ultra 350 and mineral oils | LZS 019 |
| EPDM-tubing uninsulated | 12 | 10...120 | for all bath liquids except Ultra 350 and mineral oils | RKJ 112 |
| EPDM-tubing insulated | 12 ext. dia. 35mm approx | -60...120 | for all bath liquids except Ultra 350 and mineral oils | LZS 021 |
| Silicone tubing, uninsulated | 11 | -30...100 | water, water/glycol mixture | RKJ 059 |
| Silicone tubing, insulated | 11 ext. dia. 35mm approx | -60...100 | water, water/glycol mixture | LZS 007 |
|  <ul style="list-style-type: none"> – EPDM-tube <u>not</u> for Ultra 350 and mineral oils! – Silicone oil causes pronounced swelling of Silicone rubber → never use Silicone oil with Silicone tubing! – Protect tubing with hose clips against slipping off. | | | | |

| Metal hoses single-layer insulation | Tube connection | Ø i (mm) | Ø a (mm) | Temperature range °C | Length | Ref. No. |
|---|--------------------|----------|----------|-------------------------|--------|----------------|
| MC 50 | M 16x1 | 10 | 18 | -10...400 | 50 | LZM 040 |
| MC 100 | M 16x1 | 10 | 18 | -10...400 | 100 | LZM 041 |
| MC 150 | M 16x1 | 10 | 18 | -10...400 | 150 | LZM 042 |
| MC 200 | M 16x1 | 10 | 18 | -10...400 | 200 | LZM 043 |
| Pump link | M 16x1 | 10 | 18 | -10...400 | 20 | LZM 044 |

| Metal hoses with triple insulation | Tube connection | Ø i (mm) | Ø a (mm) | Temperature range °C | Length | Ref. No. |
|--|--------------------|----------|----------|-------------------------|--------|----------------|
| MC 50 S | M 16x1 | 10 | 34 | -60...350 | 50 | LZM 046 |
| MC 100 S | M 16x1 | 10 | 34 | -60...350 | 100 | LZM 047 |
| MC 150 S | M 16x1 | 10 | 34 | -60...350 | 150 | LZM 048 |
| MC 200 S | M 16x1 | 10 | 34 | -60...350 | 200 | LZM 049 |

Further details on thermostatic liquids and hoses can be found in our special publication.

6 Unpacking, assembly and setting up

6.1 Unpacking

Goods are packed carefully, largely preventing transport damage. If unexpectedly some damage is visible on the equipment please inform the carrier or the postal authority so that it can be inspected.

Standard accessories

| | | |
|------------------------|-----------------------------|----------------------------|
| 1 Bath cover | for C 6 CP. K 6 KP. K 12 KP | Cat. No. HDQ 069 |
| 1 Bath cover | for C 12 CP. K 20 KP | Cat. No. HDQ 067 |
| 2 Bath cover | for C 20 CP | Cat. No. HDQ 067 + HDQ 068 |
| 4 Nipples 13mm dia | | Cat. No. HKO 026 |
| 4 Screw caps | | Cat. No. HKM 032 |
| 2 Closing plugs | | Cat. No. HKN 065 |
| Operating Instructions | | |

6.2 Setting up, operation as bath thermostat

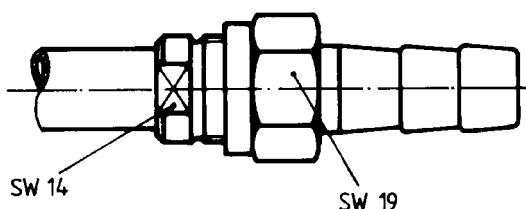
Set up the unit conveniently so that the control panel is towards the front, and ensure that the air circulation for the refrigeration system through the grille in the lower part of the unit and the ventilating openings at the back of the unit is not restricted. A minimum spacing of 20 cm between grille and wall is recommended.

Close the drain cock at the back or at the side of the bath!

When operating as bath thermostat - no external system connected up - it is advisable to ensure the internal circulation by removing the closing plugs from the pump flow and return connections. Remove the screw caps and link the pump connections together using a piece of hose. As a permanent arrangement the hose link of flexible insulated metal tubing (Cat. No. LZM 044) is the best and safest solution.

For clear-view thermostats: for improved circulation within the bath close the pump connections with the closing plugs.

Note: When loosening or tightening the screw caps (19 mm a/f), hold the threaded nipple on the tubing connections with a spanner (14 mm a/f)!



7 Connection of external systems

7.1 *Closed external circuits*

Remove the closing plugs by releasing the threaded rings (19 mm a/f) from the outflow and return connections and replace them by the tubing nipples (13 mm dia.) supplied. If the thermostat is connected to closed external circuits, additional liquid must be poured in after the thermostat has been switched on until the level in the bath remains at the correct height (approx. 2 cm below the cover plate).

At higher operating temperatures the expansion in volume of approx. 8% per 100°C during the filling must be taken into account.

For suitable tubing materials see Section 5. We recommend metal hoses for temperatures above 100°C.

With external systems at a higher level it may happen even in closed systems that the external volume drains down and the thermostat bath overflows if the pump is stopped and air enters the thermostated system!

Always ensure the maximum possible flow area in the external system (nipples, tubing, system). This results in a larger flow and therefore improved thermostatic control.

Note: Always protect the tubing with hose clips against slipping off, or use stainless steel hoses (V2A) with screwed connections.

7.2 *Open systems (baths)*

The units (except K 12 KP) are equipped with a pressure/suction pump. This can be used for the circulation of closed external circuits at higher pump outputs and in particular also of open external baths. There are two possibilities for maintaining the level in external baths:

1. The suction tubing is mounted in the external bath so that its end is at the required liquid level. The flow of the pressure stage is restricted with a tubing clamp on the hose from the pressure connection to the external bath so that flow of the pressure stage is restricted slightly below that of the suction stage. This can be noticed when air enters the suction tubing.
This operating method is not recommended, in particular at temperatures below 0°C or when using oil at elevated temperatures.
2. The preferred solution is the use of the LAUDA level controller (Cat. No. LPZ 901) which provides the functions of the adjustable level control with float, screw-on connection for external bath, and clamp fitting for 4 mm dia. Pt 100 probe.

7.2.1 Fitting the level controller

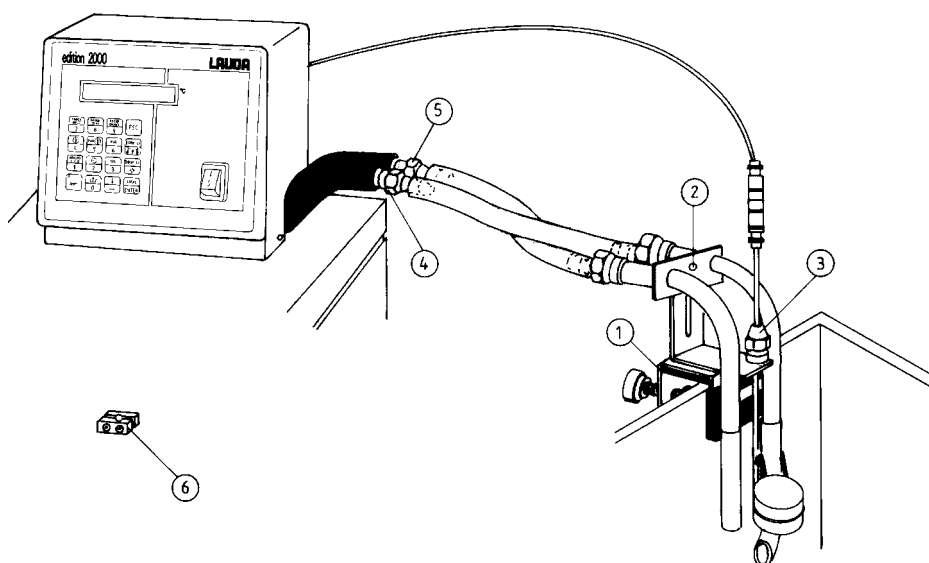
level controller

Cat. No. LPZ 901

Mount the level controller onto the external bath using the screw clamp ①. The mounting is suitable for both round and rectangular baths. Height adjustment after loosening the screw ②. Ensure a good circulation and the freedom of movement of the float!

Make the hose connections to the thermostat, connect the pressure nipple ④ to the pressure side and the suction nipple ⑤ to the suction side (see diagram).

Clamp fitting for Pt 100 probe (4 mm dia.) ③. The clamp ⑥ is not required on these units!



It is advisable to set up the external bath at the same level.

If the difference in level between the open external bath and the thermostat bath is greater than 0.5 m there is the possibility in certain applications that the control range of the level controller is not sufficient. At a higher (lower) external bath level the suction (pressure) hose should then be clamped off with a tubing clamp to such an extent that a constant level in the bath is obtained at which the float is within its control range.

Note: If the thermostat and the external bath are not at the same level it is essential to provide the venting of the connecting hoses when the pump is switched off in order to prevent overflowing.

It is preferable to use the Non-Return Fitting (see accessories) which is mounted at the highest point of the hose connection (bath or thermostat connection) and which is linked electrically to the mains output 34H.

Non-Return Fitting

Best-Nr. UD 125

Always protect the tubing with hose clips against slipping off, or use stainless steel hoses with screwed connections.

Note: When tightening the screw caps (19 mm a/f) at the tubing connections, hold the threaded nipple with a spanner (14 mm a/f)!

8 Cooling the thermostats

Because of the frictional heat of the circulating pump, thermostating without cooling can only start appreciably above the ambient temperature (see Technical data, working temperature range, lower limit). For lower temperatures it is essential to work with cooling. The following possibilities are available for cooling.

8.1 Mains water cooling

Depending on the water temperature down to 15°C. The thermostats are equipped with a cooling coil (at the rear) which is linked by tubing to the water tap and to the drain. The flow should be kept as low as possible; this saves water and improves the temperature control. Controlled cooling is possible when using a solenoid valve (see Item 3.7).

8.2 Through-flow chillers DLK 10, DLK 25 and DLK 45

They can be used, depending on the thermostat type, down to -10°C (DLK 10), -30°C (DLK 25) or -40°C (DLK 45). Use insulated hoses for the connection between the flow and return connections of the pump and the nipples of the through-flow chiller. If the thermostat operates in a closed external circuit the chiller is connected in series in the return line from the external system to the thermostat.

Always use water-glycol mixture (ratio 1:1).

9 Starting up

9.1 *Filling*

Fill the unit with bath liquid to suit the operating temperature, see Section 5. The filling volume is given under Technical data. In general the thermostat must be filled no higher than 2 cm below the cover plate. When working with thermal oils (e.g. Ultra-Therm 330 SCB) slightly less liquid should be filled in because of its expansion. The level must obviously not fall below the minimum, otherwise the low-level protection switches off the unit (see Safety circuit). The same applies to filling an external system by the pump during start-up.

Clear-View Thermostats

The unit has to be switched off and filled up to the filling marks for water or thermal oil. The mark for thermal oil, of course, only refers to the ambient temperature. When working at operating temperatures with the unit switched off the level for thermal oil can also reach the mark for water. When working with thermal oils slightly less liquid should be filled in because of its expansion. The level must obviously not fall below the minimum, otherwise the low-level protection switches off the unit (see Safety circuit). The same applies to filling an external system by the pump during start-up.

9.2 *Connection to supply*

Connect the unit only to an earthed socket (PE). Compare the details on the label with the mains voltage (see Item 4.4.2).

Model according to EMC directive EN 61326-1 (industrial areas only).*

When working without external system, ensure that the pump connections are linked together (metal hose link Cat. No. LZM 044) or closing plugs are being used.

9.3 *Basic functions*

9.3.1 *Supply switch-on*

Switch on the mains switch. The green indicating lamp lights up. The display shows consecutively

Fa. LAUDA
P-Thermostat

Type C 6 CP
V 2.XX Date

L1 Ti = 20.00°C K
L2 Ts= 10.00°C Ti I



depending on type
and software version

other values depending on
bath temperature and setpoint

* Notice only valid for EU countries!

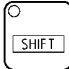


9.3.2 Standard display

Top line L1

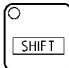


| | | |
|---|---|--|
| Ti | = | bath temperature (i = internal) |
| C | = | output in cooling range |
| H | = | output in heating range |
|  | = | cooling indication proportional to cooling actuation control |
|  | = | heating indication proportional to heating actuation control |

Bottom line L2

| | | |
|----|---|---|
| Ts | = | setpoint temperature (S = setpoint) |
| Ti | = | control variable is Ti (bath temperature), can be switched to T1 or T2 (ext. Pt 100). |
| I | = | setpoint source (I = internal = input from keys, P = programmer, R = from RS 232 C, A = analogue input socket 52 S) |

The display in line 1 (L1) can be switched by pressing the keys   and repeated operation of  to T1, T2, Ti etc.

T1. T2 = measurements of external Pt 100 probes


The display in line 2 (L2) can be switched by pressing the keys   and repeated operation of  to

Y = actual output + heating - cooling

TSI = measurement of safety comparison probe with limited resolution and accuracy


Ti. T1. T2. Ts etc..

9.3.3 Basic action on inputs and outputs

From virtually every display or input function the key  aborts and returns to the selected standard display!


Numerical inputs are always made with the SHIFT function switched off (LED in SHIFT

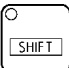
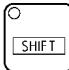

key off)!

After the last digit of a number the cursor returns to the first digit again so that corrections can easily be made before pressing the  key.

A brief beep on pressing a key means that this input is not possible!
 Error messages are indicated with text notes and accompanied by a beep. After approx. 5 sec the message disappears and the beep switches off!

9.3.4 Overtemperature switch-off point



Press the  key to indicate the current overtemperature switch-off point (To).

Important: The green LED in the  key must not be alight. If necessary bring the keypad to the basic mode by pressing the  key (LED off). If Tu is indicated, press  again to indicate To.

L2 To: _ 95.00 Ti I


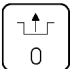
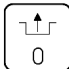
To = overtemperature switch-off point

To can be selected within the unit temperature range + 5 K.

Return to the standard display without change with . A new value is input with the number keys (SHIFT OFF); e.g. 98.70°C requires 0, 9, 8, 7, 0 . After the last digit the cursor returns to the first digit again so that any corrections can easily be made.

Pressing the key . Enter you are asked

L2 To neu Y/N 1/0__ (0)

Here the changed value must be confirmed with input  or, if no modification is requested, the initial value can be unscrambled again with input  . As default value 0 is preset.
 This interrogation was introduced in order to reach an additional protection from unintentional adjusting of the over-temperature switch-off point.

It is obviously essential to select a value above the current bath temperature (T_i) and the current setpoint (T_s), otherwise the equipment is switched off by the safety circuit with

L1 TEMPERATURE

L2 TOO HIGH!

or the message

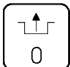
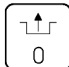
L2 T_s >> T_o

and the value is not accepted.

With activated outflow temperature limit T_{io} T_o must always be 5°C higher than T_{io} . Otherwise announcement



T_{io} >> T_o -5°C.

If the unit was switched off in the fault status an audible signal reports the stored fault when switching on.

Press the reset key  Depending on the previous sequence press  again.

If necessary check whether the overtemperature switch-off point T_o is above the current bath temperature and whether the bath is filled sufficiently!

9.3.5 Low temperature switch-off point

Press the key  (SHIFT OFF) to indicate the current low temperature switch-off point T_U . If T_o is indicated press  again to bring T_U on the display.

L2 T_U : _ -10.00 °C

T_U = low temperature switch-off point

T_U can be set up to 10 K below the working temperature range of the unit. If the bath temperature falls below T_U , T_U appears on the display in L2 so that a new value can be input if necessary. T_U operates as a setpoint limitation and as a signal.

9.3.6 Setpoint input



Press the key  (SHIFT LED OFF). L2 shows

L2 T_s :_ 20.00°C T_i I


The setpoint (T_s) can be input within the unit temperature range but not higher than the current overtemperature switch-off point.

When the input is too high, pressing  does not enter the value but instead produces the message

L2 Ts >> To


Input Ts with the number keys and the negative sign (SHIFT OFF), e.g. for -25.03°C input -, 2, 5, 0, 3 . Or for 1.93°C input 0, 0, 1, 9, 3. . After the last digit the cursor returns to the first digit again so that corrections can easily be made.

9.4 Parameter level PAR

Pressing the key  several times in the SHIFT mode (green LED in SHIFT key alight) leads successively to the input functions described below.

9.4.1 Auto-adaptation

Here it is possible to start the controller auto-adaptation by the input of 1 (SHIFT OFF)


and pressing the  key.

There should be the largest possible difference between the bath temperature and the setpoint to be entered subsequently, i.e. the time to reach the setpoint has to be longer than 5 min, preferably 10 min. In addition, auto-adaptation is obviously possible only during a heating or cooling phase which is actively influenced by the energy sources available.

Example 1: intended operating temperature approx. 70°C

1. set the setpoint to 70°C
2. within 1 minute start the auto-adaptation at the PAR level, e.g. at a bath temperature corresponding to the ambient temperature.

On reaching the setpoint the auto-adaptation switches off automatically, and the result of

the auto-adaptation can be indicated at the control parameter level  (see item 9.6).

Example 2: it is required to operate at approx. 20°C with controlled cooling.

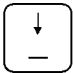




1. heat up the thermostat to approx. 60°C
2. set the setpoint to 20°C
3. start the auto-adaptation at the PAR level

9.4.2 Output limitation

Normally the maximum heating or cooling output is available. For special applications it is possible to set a limit for both heating and cooling output.

At the PAR level display select

- L1 Output
- L2 in per cent_ 100 %

Using  the display can be switched from e.g. 100%, i.e. heating output limitation, to cooling output limitation with a negative sign. SHIFT AUS. By the input of e.g. 0, 0, 5, 8  SHIFT OFF, a heating output limitation of 58% can be set. With e.g. -, 0, 9, 3  a cooling output limitation of 93% is entered. The action can be recognised by the symbols  and  flashing even at large control deviations.


Only values between 10 and 100% or -10 and -100% can be entered, otherwise the display shows the message

- L1 Output
- L2 OUT OF RANGE

9.4.3 Display resolution L1

At the PAR level display select

- L1 Display 0.001 = 1
- L2 resolution 0.01 = 0

Entering 1 (SHIFT OFF)  switches all displays in L1 to 0.001 K resolution. The temperatures are then displayed with approx. 2 digit resolution. Input "0" switches all the displays in L1 to 0.01 K resolution. Normally a resolution of 0.01 K is used.

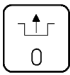

9.4.4 Contact input Fault 14 N

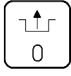
When using the contact input "FAULT" 14 N, pins 1 and 2 of the socket have to be connected together when there is no fault. If this input is not being used, a blanking plug with a link has to be plugged in. The function of the contact input fault can be switched off at the PAR level on the display

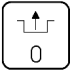
- L1 Alarm Inp. con 14 N
- L2 on = 1 off = 0

by the input of "0" (SHIFT OFF). A shorting plus is then not required.

If the alarm input has been activated in error by the input of "1", the unit can be restarted by the following inputs:

Press the key  At the PAR level select "Alarm Inp. con 14 N". Input "0" with .

Press  again.

If a fault message has been produced by opening the external signal circuit, reset by pressing the key  twice after rectifying the fault

Connections contact input "Fault" 14 N (alarm in)

3-pin flange socket to NAMUR recommendation NE 28

1 = n.o. (close)
2 = common
3 = not used

Connector plug 3-pin

Cat. No. EQS 048

Contact load approx. 5 V 2 mA. No voltage must be connected!

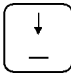

Use shielded connecting cables. Connect the shielding to the plug case. Cover the unused connectors with protective caps!

9.4.5 Baudrate RS 232

On the display at the PAR level

L1 Ser. Int RS 232

L2 Baud Rate 9600

it is possible to switch with  between 9600 and 4800. With  (SHIFT OFF) the indicated baud rate is entered.


9.4.6 Menu language

On the display at the PAR level

L1 Lang. Germ = 0

L2 Engl. = 1 French. = 2

the menu language can be selected. Enter the corresponding code numbers 0, 1 or 2 with

 (SHIFT OFF).

9.4.7 Calibrating the analogue output channels

The 90% values of the analogue voltage outputs channel 1 and 2 or the analogue current output of channel 1 can be calibrated separately for channel 1 (voltage or current) and channel 2 (voltage). The factory calibration on channels 1 and 2 for 0...10 V = - 100...400°C is performed at 9 V = 350°C.

In special cases, e.g. to correct scaling deviations of instruments connected to the output, or if channel 1 is to be a current output, the output can be calibrated by the user.

At the PAR level display select

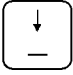
L1 Analogue outp. Cal ?


L2 Chan 1 = 1 Chan 2 = 2


Input SHIFT OFF 1 , or 2 for channel 2..

Depending on the selected configuration of the analogue outputs (see Item 9.12) the socket 52 S (analogue signals, see Item 9.10) at pin 2 carries a voltage signal of approx. 95% or 9.5 V, or pin 5 the corresponding current signal of approx. 19 mA in case of current configuration for channel 1.

Using a precision multimeter or e.g. a temperature recorder set the output signal to 9 V or

18 mA or the corresponding temperature by the repeated operation of the key  (SHIFT ON).

Pressing  leaves the menu and the most recent value setting is entered.

If the value was selected too low, leave the PAR level with  and make a new selection.

The calibration of channel 2 is similar. Connect the measuring instrument to pin 1 (voltage signal only).

9.4.8 Operation with through-flow chiller DLK 45 with proportional cooling

At the PAR level display select

L1 DLK normal = 0

L2 DLK auto = 1


Select „DLK normal“ by pressing 0 in order to have the function as already known of the two outputs 19 H and 34 H.

Select „DLK auto“ by pressing 1 in order to operate a through-flow chiller DLK 45 with automatic compressor control and proportional cooling; see also operating instructions of the through-flow chillers.


9.5 Calibration of the temperature measurement circuits

With the calibration function the indications of the three temperature measuring points bath temperature T_i , external Pt 100 probe T_1 and external Pt 100 probe T_2 can be set to a known accurate value. The resulting correction is processed additively over the entire temperature range.

Check first that a sufficiently accurate reference is available, otherwise it is better to use the factory calibration which gets lost by overwriting!

Pressing the key  in the SHIFT mode (green LED in shift key alight) produces the display

L1 CALIBRATE
 L2 $T_i = 0 \quad T_1 = 1 \quad T_2 = 2$

The channel to be calibrated is selected with 0, 1 or 2 .

When selecting an unused channel, e.g. if Pt 100 on T_2 is not connected, the display shows


L1 Ext Pt 100 not
 L2 connected

For calibration a sufficiently accurate reference temperature measurement should be possible, and the measurement point temperature should be constant.

The display shows

L1 $T_1 \quad 61.04^\circ\text{C}$
 L2 $T_c _ . \quad ^\circ\text{C}$

The value shown in L1 is the measured value obtained without any correction using probe and electronics without calibration. Now enter the real value for the measurement point T_1 (e.g. 60.00°C).

Example: 0. 6. 0. 0. 0 

T_i or T_2 can be calibrated in the same way.

In order to avoid dangerous conditions the correction is limited to $\pm 5 \text{ K}$. In case of larger corrections the display shows


L1 CORRECTION VALUE
 L2 TOO LARGE

and the entered value is not accepted.

You can leave the calibration level with .

9.6 Control parameters

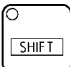
9.6.1 Indication and input of the control parameters

Pressing the key  several times in the SHIFT mode (green LED in shift key alight) shows the outflow temperature limitation, the correction limitation and the control parameters X_p, T_n and T_v on the display in L2.


Example:

| | | | | |
|----|-------------------|---------------|----|---|
| L2 | T _{io} : | <u>1</u> 20°C | Ti | I |
| L2 | T _d : | 30°C | Ti | I |
| L2 | X _p : | 0.5°C | Ti | I |
| L2 | T _n : | 12.0 s | Ti | I |
| L2 | T _v : | 2.0 s | Ti | I |

In order to use control parameters other than those found by auto-adaptation (see Item 9.4.1) the values can be entered in the appropriate display after switching off the

SHIFT function, pressing . LED off.

Example for X_P:

0. 0. 1. 0. 

if the required value is 1.0°C.

For values above 200.0°C or 200.0 sec the message

L2 OUT OF RANGE

Appears.

9.6.2 Recommendations for the control parameters

In most cases satisfactory control results are obtained with the following control parameters:

| bath liquid | water | oil |
|----------------|-------|------|
| X _P | 0.5°C | 1°C |
| T _n | 10 s | 25 s |
| T _v | 2 s | 5 s |

9.6.3 Bath temperature limitation


The limitation of the bath temperature is an additional warning and switch-off function switching off the heating at a selectable value; i.e. the heating output is set to "0". This protects the unit from a continuous cut-off via the safety circuit especially during external control at certain operating conditions.

To enter the switch-off point T_{iO} proceed as described in Item 9.6.1, and switch the display to input and indication.

Example:

L2 T_{iO} 120 °C

Change the value by entering numerical inputs with a resolution of 1°C. Enter the new

value by pressing .

It is possible to select values within a range from 50°C to the selected overtemperature switch-off point T_O -5°C. If this range has not been respected the display shows the message

L2 $T_{iO} >> T_O - 5$ °C

Of course T_{iO} has to be set above the setpoint T_S ; otherwise the display shows the message

L2 $T_S > T_{iO}$

The bath temperature limitation can be switched off by entering

L2 T_{iO} 000

If the bath temperature T_i exceeds the selected switch-off point the display shows

Example:

L2 T_{iO} 120 °C

and there is an acoustic signal. The heater switches off. As soon as the temperature has dropped the unit starts working again.

9.6.4 Correction limitation

During the operation with external control it may be necessary not to exceed the difference between the bath temperature T_i and the measuring point for the external control T_1 or T_2 , e.g. in order to get a smooth heating of the material or the vessel.


Such a value can be selected by the variable T_d . If the value T_d is exceeded the heating or cooling output is set to "0". If this function is activated the times for heating up or cooling down may be extended.

To enter the difference value T_d , proceed as described in Item 9.6.1 and switch the display to input and indication.

Example:

L2 T_d _ 30°C

Change the value by entering numerical inputs with a resolution of 1°C.

Enter the new value by pressing .

It is possible to select values within the temperature range from 5°C to 150°C. If this range has not been respected the display shows

L 2 OUT OF RANGE

and there is an acoustic signal.

This function can be switched off by entering

L2 T_d 000 °C

9.7 External control

9.7.1 External measurement inputs and external controller

The units have two Pt 100 temperature measurement inputs whose measurements can be indicated (T_1 , T_2).

You can connect the external Pt 100 (T_1 , T_2) at the rear connectors 10 S in 4-wire circuit.

Use shielded connecting cables. Connect the shielding to the plug case. Cover the unused connectors with protective caps!

Pin connections sockets 10 S Pt 100

pin

| | | | |
|---|---|---|--------------|
| 1 | + | I | current path |
| 2 | + | U | voltage path |
| 3 | - | U | voltage path |
| 4 | - | I | current path |



Pt 100
DIN IEC 751


Plug, 4-pin Lemos, for Pt 100 connection


Cat. No. EQS 022

One probe can be selected for the actual value for external control. The unit then operates with cascade control to this actual value, i.e. the unit controls the temperature at the external measurement point to the selected setpoint by suitably altering the bath temperature. Thus the influence of disturbances (changes of load or through-flow, etc.) can be reduced considerably or eliminated totally.

9.7.2 Start of external control

Connect platinum resistance thermometers to both of the Pt 100 inputs 10 S (T₁ and T₂). It is sensible but not essential to use T₁ if only one input is in use.

The external control with the measuring point T₁ is switched on with the key  in the SHIFT mode; pressing the key again switches to T₂ as control variable.

Pressing the key  once more (SHIFT ON) switches back to the bath control (internal) T_i.

In L2 the position before the final one shows the parameter used as control variable.

```

L2      Ts      = 120.35   Ti   I
                                T1
                                T2
  
```

This setting remains stored in case of a fault or after the power is switched off.


If T₁ is selected but no probe has been connected the message

```

L1      Ext Pt 100 not
L2      connected
  
```

appears.

The unit then switches the control variable to T₂ automatically. If T₂ is also not connected, the thermostat switches to T_i.

When changing the setpoint for more than 10°C it may be possible to achieve an improved control result by restarting the external control from the control variable T_i (with SHIFT ).

After power OFF the unit operates with control from the bath (T_i) for safety reasons; external control must be re-selected after power ON as explained above.

9.7.3 Notes

When operating with external control it is essential to ensure that the probe for the control variable is in good thermal contact with the liquid, otherwise a poor control result must be expected, or the control may be completely ineffective.

Proceeding from the control parameters used for bath control the control may have to be adapted either by auto-adaptation (see Item 9.4.1) or by the input of the control parameters.


Important: set the overtemperature switch-off point T_O (see Item 9.3.4) sufficiently high since the bath temperature may under certain circumstances become much higher than the setpoint.

9.8 **Working with controlled cooling**

Operation with controlled cooling requires a solenoid valve.

Insert the plug of the solenoid valve into the socket (19 H) on the back. The solenoid valve can be fitted either on the cooling coil or on a 1/2" water tap. Although fitting on the cooling coil is the usual method, mounting directly on the water tap is preferable for two reasons: When the valve is closed, the connection hose to the cooling coil is not under pressure; therefore there will not be a pressure surge when the valve is switched on, and the danger of the hose bursting is much reduced. Use hose clips!

With controlled cooling operation the solenoid valve switches with a cycle time of about 6 s. L1 indicates on the right the symbol for the status of the solenoid valve. Restrict the water flow as much as possible at the water tap. This produces improved control and saves cooling water.

Note: Ensure that the cooling coil connections  are used. Do not mix them up with the pump connectors!

It is essential to ensure free outflow from the cooling coil, especially at operating temperatures above 100°C because of steam formation! The use of controlled cooling is particularly helpful when initiating exothermal reactions or in programmer operation.

Solenoid valve for cooling water control

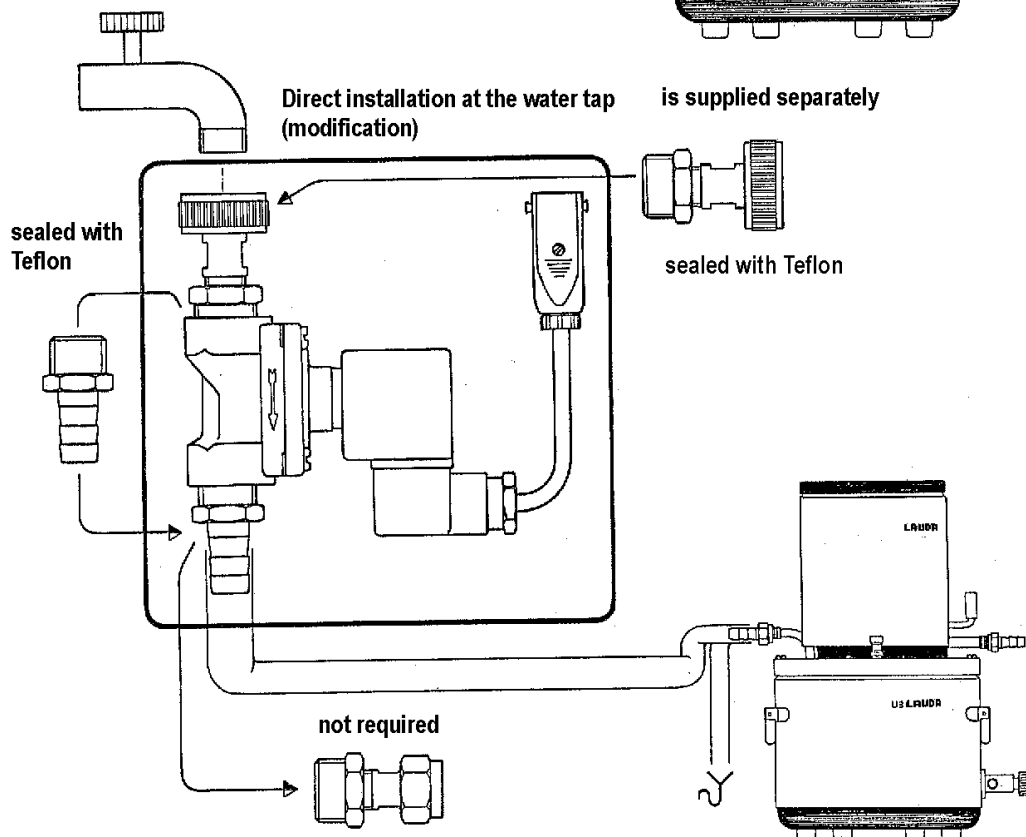
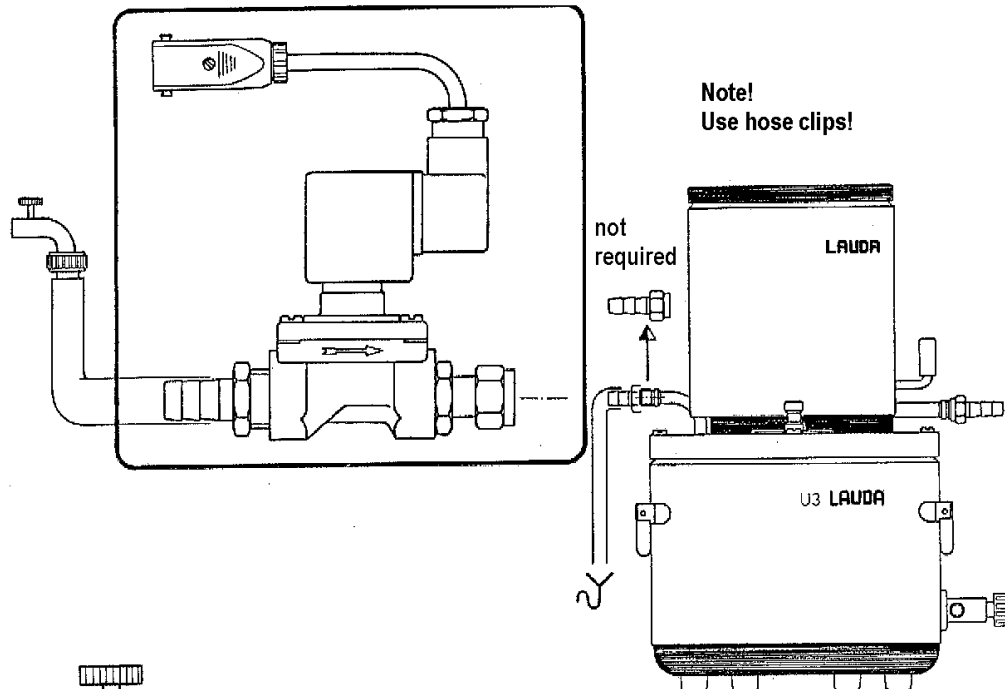
Cat. No. UD 085

Mating plug for other solenoid valve

Cat. No. EQS 005

Mounting instructions for solenoid valve

Direct installation on the cooling coil (as supplied)



9.9 Operation with programmer

Temperature programmes with up to 99 segments can be stored and processed. A segment consists of a target temperature which is to be reached at the end of the segment, and the time duration of the segment. The time "00:00" for temperature differences is possible. In connection with the tolerance range monitoring the programme continues not until the target temperature is reached. It is useful to prepare a time-temperature diagram before programming and to check whether the energy balance enables the programme speed.


Set the overtemperature switch-off point To to a value slightly above the highest bath temperature to be expected (see Item 9.3.4).

9.9.1 Programme input

Press the key  in SHIFT mode. The display shows

L1 PROG. INP
 L2 Tstart: _ . °C

Enter here the starting temperature of the programme. SHIFT OFF (automatically), e.g.

for 60.00°C input 0,6,0,0,0 

The display shows

L1 PROG. INP SEG.01
 L2 T : _ . °C : h


Now enter the target temperature and the time for the first segment, e.g. for 140.00°C in

the time 2 h 00 min 1, 4, 0, 0, 0  then 0,2  then 0, 0 

The display shows

L1 PROG. INP SEG.02
 L2 T : _ . °C : h

Now enter the target temperature and the time for the second segment, e.g. for a phase at a constant temperature 140.00°C and 1 h 30 min.

After the last programme segment press the key  once more.

The display shows

L1 PROG. INP
L2 NO OF CYCLES:_

Input 1 ... 99 is possible.

With more than one cycle it is convenient to have the final temperature and the starting temperature T_{start} at the same level!

Afterwards a tolerance range can be input for monitoring the programme.

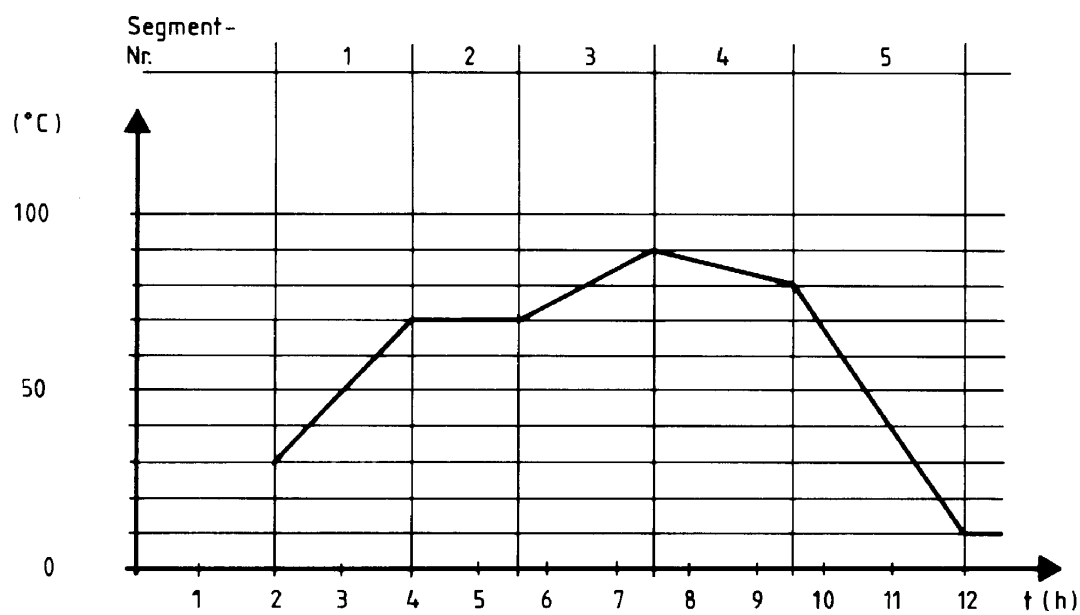
The display shows

L1 PROG. INP
L2 TOL. RANGE _.

Now you can input a tolerance range value from 0.1 to 9.9°C. I.e. if the control variable (bath temperature or external temperature T_1 or T_2) deviates from the set temperature of the segment by more than the tolerance range value while the programme is running, the programme sequence will be stopped until the control variable is within the tolerance range again. At the same time a "T" appears on the right in L2


The input of 0.0 switches off the tolerance range function.

9.9.2 Example of a programme




| Segment-No. | | Input | Press <div>LOCAL ENTER</div> |
|-----------------|-----|---------------------|----------------------------------|
| Tstart | | 03000 | 1x |
| 1 | °C | 7000 | 1x |
| | h | 02 | 1x |
| | min | 00 | 1x |
| | | | |
| 2 | °C | 7000 | 1x |
| | h | 01 | 1x |
| | min | 30 | 1x |
| | | | |
| 3 | °C | 9000 | 1x |
| | h | 02 | 1x |
| | min | 00 | 1x |
| | | | |
| 4 | °C | 8000 | 1x |
| | h | 02 | 1x |
| | min | 00 | 1x |
| | | | |
| 5 | °C | 1000 | 1x |
| | h | 02 | 1x |
| | min | 30 | 2x |
| | | | |
| Cycles | | 1...99 | 1x |
| Tolerance range | | ±°C (0.0) 0.1...9.9 | 1x |

9.9.3 Programme test

After the input of the programme it is advisable to check that the programme buffer contains the correct data. This is done with the key  in the SHIFT mode.


Pressing the key repeatedly produces the same sequence as during the input of the programme.

9.9.4 Changing the programme data

Select the data line to be changed as in "programme test", SHIFT OFF. This resets the data of the indicated segment. Then the data can be input as usual. Enter the new data with  each.





9.9.5 Programme start, interruption and abort

It is useful to bring the operating temperature of the thermostat to the programme start temperature T_{start} before starting the programme, or to automatize it by the tolerance range function.

Start the programme sequence with  in the SHIFT mode. The programme sequence can then be followed by the indication of the setpoint T_s .

L1 shows on the right the segment number, and L2 shows as a setpoint source a P for "programme" on the right.

The keys 1 to 9 are blocked while the programme is running.

The sequence of the programme can be stopped with SHIFT  and then be re-started with SHIFT  During the interruption of the programme the display shows a W (wait) on the right in L2. The programme can be aborted with  and then SHIFT  within 2 sec. Afterwards the programme can only be started with segment 1.

9.10 Connection for analogue signals socket 52 S

6-pin flange socket according to NAMUR recommendation NE 28.

| | |
|--------|--|
| Pin 1: | voltage output temperature signal channel 2: setpoint T_s , bath temperature T_i , external Pt 100 T_1 or T_2 can be selected. Scaling can be as follows: 0...10 V corresponding to a temperature range selected within the working temperature range (e.g. 50...80°C) minimum load 4 kOhm. or 0...6 V = -200...400°C = 10 mV/K 0°C = 2 V |
|--------|--|

| | |
|--------|--|
| | or 0... 10 V = -100...400°C or 0...10 V = 0...100°C |
| Pin 2: | voltage output temperature signal channel 1, other data as pin 1 |
| Pin 3: | ground for all signals |
| Pin 4: | setpoint voltage input; scaling can be selected as pin 1. R _i = 12 kOhm approx. (+ pin 4; - pin 3) |
| Pin 5: | current output temperature signal channel 1; signal selection as pin 1. Can be configured for 0...20 mA or 4...20 mA. Scaling can be: 0...20/4...20 mA = -100...400°C or 0...20/4...20 mA = 0...100°C or 0...20/4...20 mA = a temperature range selected within the working temperature range (e.g. 50...80°C) maximum burden 330 Ohm. |

Connect only either pin 2 or pin 5!!

| | |
|--------|---|
| Pin 6: | setpoint current input; configuration and scaling as pin 5. Burden 320 Ohm approx. Maximum voltage 15 V! |
|--------|---|

Connector plug, 6-pin

Cat. No. EQS 057

Use shielded connecting cables. Connect the shielding to the plug case. The mass for all signals (pin 3) must not be connected with ground! If a connection to the ground cannot be avoided use a potential-free signal bridge in between.
 Cover the unused connectors with protective caps!

9.11 *Analogue inputs*

A setpoint in the form of an analogue current or voltage signal can be provided by connection to the socket "Temp.-Signal" 52 S (see Item 9.10).

- By pressing the key  in the SHIFT mode the display shows

L1 ANALOG INP/OUTP


L2 INP = 0 OFF = 1_

Pressing 0  selects the configuration and scaling of a setpoint input.


The display shows

L1 ANALOG INP. OFF = 0


L2 ON = 1 CONF = 2_

- Input 1  switches in a previously configured input as setpoint, and L2 shows an "A" at the right end, indicating that the setpoint is determined by the analogue input.

This condition remains stored in case of a fault or after the power is switched off.

Input 0  switches the setpoint back to setpoint source I internal, i.e. key input.

Scaling takes place interactively by applying the voltage and current values corresponding to the appropriate temperature range limits to the corresponding input. For pin connections for voltage or current input on socket 52 S see Item 9.10. This method compensates various scaling errors, e.g. also those of the sources connected.

- Pressing 2  configures and scales the setpoint input. The display shows

L1 ANALOG INPUT

L2 U = 0 I = 1

Select a voltage range with 0  Voltages in the range 0...10.5 V can be handled.

A current range is selected with 1  Currents in the range 0...22 mA can be handled.


- The display shows

L1 ANALOG INPUT

L2 T_{min} = _ . °C

Input the lowest temperature of the range which corresponds to the lowest voltage or current value of the range to be scaled.

Example: range 0...120°C should correspond to 0...10 V approx.

Input 0, 0, 0, 0, 0 .

The display shows

L1 ANALOG INPUT

L2 $T_{\max} = _ . \text{ } ^\circ\text{C}$



Input the upper limit of the temperature range 1, 2, 0, 0, 0 .

If the current input is selected, the programme asks whether 0...20 mA or 4...20 mA is required.

The display shows

L1 CURRENT INPUT

L2 0-20 = 0 4-20 = 1


Select 0  or 1  This menu item is omitted when the voltage input has been selected.

The display shows

L1 INPUT CAL.?

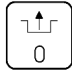
L2 YES = 1 NO = 0

Here the decision is made whether an automatic calibration procedure is started, or whether the voltage or current values from the last calibration procedure are retained with

input 0 .

The display returns to the standard display. With a new calibration the voltage or current source (e.g. setpoint unit, programmer) must be connected up. The range limits must be adjustable.

The unit may switch to fault if the input signal is not connected. If this is the case it is


necessary to connect the input signal first. Then press the reset key  and calibrate in the same way as described above.

Select recalibration with 1 .

The display shows

L1 SET U_{min}

L2 YES = 1

When the voltage or current corresponding to the lower range limit is applied to the input, confirm this by input of 1 .

The display shows


L1 -----wait-----

L2 XXXXXXXXXXXXX

The calibration takes approx. 20 sec. Then the display shows

L1 SET U_{max}

L2 YES = 1

When the voltage or current corresponding to the upper range limit is applied to the input, confirm this by input of 1 .

The display shows

L1 -----wait-----

L2 XXXXXXXXXXXXX

The calibration takes approx. 60 sec. Then the display returns to the standard display.
The calibration is finished.


- Switch-on the external setpoint from an analogue input as described in Item 9.11.2.

9.12 Analogue outputs

Two analogue output channels are available at the socket "Temp.-Signal" 52 S (see Item 9.10). They can be set to carry the temperature values

| | | |
|----------------|---|---|
| T _i | = | bath temperature |
| T ₁ | = | temperature at ext. Pt 100 T ₁ |
| T ₂ | = | temperature at ext. Pt 100 T ₂ |
| T _s | = | setpoint |

9.12.1 Temperature signal channel 1


Channel 1 can be configured at socket 52 S on pin 2 as voltage output or on pin 5 as current output. Press key  in SHIFT mode.

The display shows

L1 ANALOG INP/OUTP
L2 INP = 0 OUTP = 1

Select the processing of the outputs by pressing 1 . The display shows


L1 Analog outputs
L2 Chan 1=1 Chan 2=2

Select channel 1 by pressing 1 . The display shows

L1 Analog output
L2 U = 0 I = 1_

Select the current output with 1 . The display shows

L1 CURRENT OUTPUT
L2 0-20=0 4-20=1

Select the required current range 0...20 mA or 4...20 mA by pressing 0 or 1 

The current range selection is omitted if the voltage range 0...10 V has been selected in the previous menu. The display shows the scalings available for selection

L1 CONFIGURABLE.=1

L2 analog output _



By pressing the key (SHIFT ON) the pre-set scalings are displayed consecutively.

The selection is made by the input of the appropriate code (SHIFT OFF).

Configurable means that the temperature range required to correspond to the voltage range 0...10 V, current range 0...20 mA or 4...20 mA can be determined by setting the range start (T_{min}) and the range end (T_{max}).

The display shows


L1 CONFIGURABLE.=1

L2 $T_{min} = _ . \text{ } ^\circ\text{C}$

Example: Bereich 20...220°C

Input 0, 2, 0, 0, 0 . The display shows

L2 $T_{max} = _ . \text{ } ^\circ\text{C}$

Input 2, 2, 0, 0, 0 .

The following fixed scalings are available:

-200...400°C = 0...6 V = 10 mV/K


0°C = 2 V Code 2


-100...400°C = 0...10 V oder 0...20 mA or 4...20 mA Code 3

0...100°C = 0...10 V oder 0...20 mA or 4...20 mA Code 4

The display then shows

L2 $T_i T_1 T_2 T_s 0-3_$

E.g. to set the bath temperature on channel 1 input 0 .

Similarly for the temperature signal of the external Pt 100 T_1 input 1  etc. The display then returns to the standard display.

9.12.2 Temperature signal channel 2

Channel 2 is purely a voltage output at socket 52 S on pin 1. The selection is made as described for channel 1 in Item 9.12.1 except that current ranges cannot be selected.

9.13 Safety function

The operation of the safety devices of the units has already been described under Item 4.1.

After starting up the user should confirm the correct operation of the safety devices. If the unit operates without supervision we recommend that this check should be carried out daily.

9.13.1 Low-level cut-out

For a correct operation of the low-level cut-out it is essential that the float switch operates correctly. To check this, lower the level in the bath by draining away some of the liquid. When the bath level falls below the minimum level (approx. 20 mm above the upper heater winding) the pump, the heating and the refrigeration unit switch off on all poles.

The display shows the message

L1 LEVEL
L2 TOO LOW

and there is a warning beep.

To restart fill up the bath and press the reset key  twice (with approx. 1 sec interval).

9.13.2 Adjustable overtemperature limiter

To check it the switch-off point TO has to be set below the current bath temperature. Note that an input of TO below the setpoint TS produces the message.

L2 TS >> TO

and the previous value for TO is retained. Therefore the setpoint TS usually has to be lowered first by a few degrees before carrying out this test. The overtemperature switch-off point can then be set e.g. 1 K below the current bath temperature.


Example:


| | | |
|----|---|------|
| Ti | = | 60°C |
| TS | = | 60°C |
| TO | = | 65°C |

To check the operation of the overtemperature limiter, press  (SHIFT OFF). Input

Ts = 20°C, press , the display shows

| | | | |
|----|----|---|--------------|
| L1 | Ti | = | 60.00°C C |
| L2 | To | : | 65.00°C Ti I |

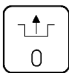
If Tu is shown in L2, press  again!


Now input 0, 5, 8, 0, 0 .


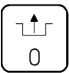
The switch-off point of the overtemperature limiter is now 2 K below the bath temperature Ti. The display shows the message

| | |
|----|-------------|
| L1 | TEMPERATURE |
| L2 | TOO HIGH! |

with a beep. Heating and pump are switched off on all poles.

To restart the unit press the reset key . The display shows the standard display.

Now press  and set TO to a value above the bath temperature, e.g. 70°C:

Input 0, 7, 0, 0, 0 . Then press the key  once more. The unit returns to normal operation.

Note: The overtemperature switch-off point has to be set at least 25 K below the fire point of the bath liquid used according to EN 61010.

In case of any malfunction in Items 9.13.1 and 9.13.2 the unit has to be taken out of use immediately and checked by an engineer, otherwise its safety is no longer ensured.

9.13.3 Connection potential-free contact "Combined fault" 12 N (Alarm off)

3-pin flange connector conforms to NAMUR recommendation NE 28

1 = n.o. (make)

2 = common

3 = n.c. (break)

1,2 are linked when unit operation is OK

Use shielded connecting cables. Connect the shielding to the plug case. Cover the unused connectors with protective caps!

Coupling socket 3-pin

Cat.No. EQD 047

9.14 Serial interface RS 232 C

-pin sub-D socket 53 S

9.14.1 Data of the RS 232 C interface

Cables used (computer end)

| <u>Computer</u> | | <u>Thermostat</u> | | |
|-----------------|-------|-------------------|-------|---------------------------|
| 25-pin | | 9-pin | 9-pin | |
| 3 | R x D | 2 | 2 | T x D (Transmitted Data) |
| 2 | T x D | 3 | 3 | R x D (Received Data) |
| 7 | SG | 5 | 5 | Signal Ground |
| 6 | DSR | 6 | 6 | DTR (Data Terminal Ready) |
| 4 | RTS | 7 | 7 | CTS (Clear to send) |
| 5 | CTS | 8 | 8 | RTS (Request to send) |

Use shielded connecting cables. Connect the shielding to the plug case. Cover the unused connectors with protective caps!

Using this interface it is possible to transfer the following data from or to a computer with a suitable interface:

1. Transfer of the setpoint from the computer to the thermostat
2. Read-out of the bath temperature T_i , the external temperature T_1 , the external temperature T_2 and the setpoint on the unit
3. Transfer of low temperature and overtemperature switch-off point
4. Read-out of the set overtemperature and low temperature switch-off point
5. Read-out of the fault signal
6. Transfer of the ramp segments and their processing
7. Status signal

8. Read-out of the control parameter and transfer
9. External controller status and start



9.14.2 General principles



The interface operates with two stop bits, no parity bit and with 8 data bits. The transfer rate can be set to 4800 baud or 9600 baud (see Item 9.4.6).

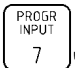

Values from the computer can be transferred directly to the thermostat, i.e. transmitted, e.g. OUT, SEG and START commands, or data can be transmitted from the thermostat to the computer on request with an IN command. An OUT, SEG or START command, if transmitted correctly, is always acknowledged by the thermostat with the message "OK" followed by LF and CR.

This message, like any other response, has to be requested by the computer!

Any output instruction (OUT, SEG, START, STOP) switches the thermostat to remote operation. This can be recognized by an R (setpoint source RS 232) on the right in L2.

Then all the keys are locked except for the functions "SHIFT  " and "SHIFT  ". If there is no output instruction from the connected computer the keyboard can be activated

until the next output instruction by pressing the keys  . The data requests by

the thermostat (IN commands) only lock the programme keys "SHIFT  " and "SHIFT  ". All the other key functions are in operation. In the following text the symbol "_" will be taken to mean blank (no character).

RS 232 interface and controller are operated by a single processor; for optimum control it is therefore advisable to have pauses of at least 100 msec between the interface commands.

9.14.3 Output commands

OUT_XXX.XX

Setpoint transfer with up to 3 places before the decimal point and up to 3 places behind. This includes the negative sign. Transfer can take various forms, e.g. for 5.00°C: 005.00, 05, 05.0, 005, 5.00.

A BASIC programme for the IBM PC which can be used to transfer any values between the set upper limit (see Item 9.3.4) and which displays the response "O.K." or a possible error message, may be as follows:

Note: set baud rate to 4800 (see Item 9.4.5)!

```

10 OPEN "COM1:4800,N,8,2" AS #1
20 CLS
30 LOCATE 8,5:PRINT SPC(70)
40 LOCATE 8,5
50 INPUT "Enter your command (without OUT_);VALUE$
60 PRINT #1;"OUT_"+VALUE$
70 INPUT #1;A$
80 LOCATE 12,5:PRINT SPC(50)
90 LOCATE 12,5:PRINT "Response of the thermostat";A$
100 TI = TIMER+1
110 IF TI > TIMER THEN 110
120 GOTO 30
130 END

```

The following values can be transmitted similarly to the thermostat:

| | |
|---------------------|---|
| OUT_LXXX.XX | switching point for low temperature (usually set to the lower range limit of the thermostat) |
| OUT_HXXX.XX | overtemperature switch-off point. For safety reasons it is essential that, after the transfer, this value is read back with the command IN_9 and checked! |
| OUT_XPXXX.XX | setting of the control parameter X_p for the controller |
| OUT_TNXXX.XX | setting of the control parameter T_n |
| OUT_TVXXX.XX | setting of the control parameter T_v |
| OUT_RT1 | switches the control variable to the source external Pt 100 T ₁ (external control) |
| OUT_RT2 | switches the control variable to the source external Pt 100 T ₂ (external control) |
| OUT_RT _i | switches the control variable to the source T _i (probe in the bath); control according to the bath temperature |
| SEG_XXX.XX_XX:XX | using this programme segment command a segment can be written into the programmer buffer. It indicates the target temperature and the segment time hours (2 digits max.) and minutes (59 max.). The segment start is formed by the current setpoint, i.e. before the transfer of a programme segment it is useful to transfer a setpoint as a segment start suitable for the subsequent segment, using OUT_XXX.XX. "_ " blank (no character). |

| | |
|--------------------|---|
| SEG_(XX)_XXX_XX:XX | Single segment with segment number, used when whole temperature programmes are to be loaded from the computer to the thermostat. Thus, in contrast to the command SEG_, several segments may be transmitted. The programme starts with the latest setpoint; therefore please check before START whether the setpoint suitable for the first segment is available in the unit. |
| OUT_TBX.X | the tolerance range value is $\frac{1}{2}$ of the value of the total range; i.e. 0.5 is 0.5 K; range 0.1...9.9 K. 0.0 switches off the tolerance range function. |
| OUT_CYXX | number of the programme cycles, range 1...99. 0 switches off the function, i.e. the programme is repeated until it is stopped manually. |
| START | starts the segment contained in the programme buffer |
| STOP | stops the programme segment run. With START the programme segment starts again from the beginning. |

9.14.4 Requesting data from the thermostat

| | |
|-----------------------|---|
| IN_1 | indication of the bath temperature (Ti), i.e. the request of the thermostat to transmit the bath temperature. |
| IN_2 | indication of the temperature value at the external probe T1 |
| IN_3 | indication of the current setpoint (Ts) |
| IN_4 | status signal, 7 characters |
| char 1 from the left: | overtemperature fault = 1, no fault = 0 |
| char 2: | low level fault = 1, level OK = 0 |
| char 3: | programmer segment running = 1, programmer segment off = 0 |
| char 4: | control according to the bath temperature (Ti) = 0, T1 = 1, T2 = 2 |
| char 5: | setpoint set by analogue inputs = 1, analogue inputs off = 0 |
| char 6: | indicates whether external Pt 100 T1 is connected = 1, or not connected = 0 |
| char 7: | indicates whether external Pt 100 T2 is connected = 1, or not connected = 0 |

IN_5 invalid

IN_6 invalid

IN IN_7 indication of the temperature of the external probe T2

IN_8 indication of the current low temperature switch-off point TU

IN_9 indication of the current overtemperature switch-off point TO

IN_A indication of the current value of Xp

IN_B indication of the current value of Tn

IN_C indication of the current value of Tv

A BASIC programme used to transfer values from the thermostat to the computer and to display them specifying the channel number (e.g. 1 for IN_1, i.e. the bath temperature), is as follows:

Note: set baud rate to 4800 (see Item 9.4.5)!

```

10  OPEN "COM1:4800,N,8,2" AS #1
20  CLS
30  LOCATE 8,5:PRINT SPC(20)
40  LOCATE 8,5
50  INPUT "Channel No.";NO$
60  PRINT #1;"IN_"+NO$
70  INPUT #1;A$
80  LOCATE 12,5:PRINT SPC(50)
90  LOCATE 12,5:PRINT "Response of the thermostat";A$
100 TI = TIMER+1
110 IF TI > TIMER THEN 110
120 GOTO 30
130 END
    
```

The isolation of the status data may be as follows:

| | | |
|------------------|---|------------------------------|
| LEFT \$ (A\$,1) | = | overtemperature fault |
| MID \$ (A\$,2,1) | = | low-level fault |
| MID \$ (A\$,3,1) | = | programme segment running |
| MID \$ (A\$,4,1) | = | control by Ti, T1 or T2 |
| MID \$ (A\$,5,1) | = | analogue input on/off |
| MID \$ (A\$,6,1) | = | external Pt 100 T1 connected |
| RIGHT \$ (A\$,1) | = | external Pt 100 T2 connected |

9.14.5 Error messages on the computer

The following error messages can be reported from the thermostat to the computer during operation:

- ERR-2: invalid inputs (e.g.: overflow of the input buffer)
ERR-3: invalid command
ERR-5: invalid command when switching the control variable for the controller, e.g. external controller OUT_RT2. Other command than OUT_RTI, OUT_RT1, OUT_RT2.
ERR-6: temperature value cannot be set
ERR-7: syntax error in channel number
ERR-8: channel does not exist

9.15 Timing clock function

The unit is equipped with a clock indicating day, month, year, weekday, hours and minutes and provides these for the timing function.

The clock is backed for approx. 10 years by a built-in battery so that the clock continues to operate even when the unit is not connected to the electrical supply.


9.15.1 Setting and indication of date and time


This function is only required when changing from summer to winter time and vice versa, or when the unit is being operated in other time zones. Date and time are set at the factory when the unit is started up for the first time.

By pressing the key  in the SHIFT mode the display shows

L1 clock = 0 activ = 1

L2 SET = 2 FUNCT = 3

By the input of 0  date and time are indicated.  returns to the standard display.

The input of 2  (SET) allows date and time to be altered. The display shows:

L1 DA MO YE H MI

L2 _ . . :


Day, month, year, weekday, hours and minutes are now input in sequence.

Weekday code.

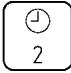

- | | |
|---|-----------|
| 1 | Monday |
| 2 | Tuesday |
| 3 | Wednesday |
| 4 | Thursday |
| 5 | Friday |
| 6 | Saturday |
| 7 | Sunday |

The hours are input from 0 to 24 (factory-set to Central European Time).

Example: L2 19.01.94 3 16:05

Terminate the input with .


9.15.2 Timing clock function

By pressing the key  in SHIFT mode the display shows the menu as described in Item 9.15.1. Select the timing function FUNCT = 3 with 3 .

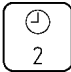

The display shows

- | | |
|----|-------------------|
| L1 | Thermostat ON = 1 |
| L2 | OR OFF = 0 |


Here it can be selected whether the thermostat should automatically switch on or off at the time to be selected subsequently in Item 9.15.3. Input either 1 or 0 as appropriate. A display to input date and time appears.

Input here the switching point for the timing function as described in Item 9.15.1 and enter it with . The display returns to the standard display.

9.15.3 Activating the timing clock function

By pressing the key  in SHIFT mode the display shows the selection menu as described in Item 9.15.1. Activate the timing function with ACTIV = 1 by the input of 1 . Again the display shows



- | | |
|----|----------------|
| L1 | Clock |
| L2 | ON = 1 OFF = 0 |

Normally 1  is input here to activate the timing function.

The timing clock symbol now appears in line L1 before the last position. If the previous selection was that the thermostat should switch on automatically, the thermostat now switches off and starts up at the selected time.

When the unit has been switched off through the timing function, the display shows

L2 CLOCK STOP!

The activated timing function can always be switched off with off = 0, i.e. with 0 . In addition the timing function can be cancelled at any time with .

10 Maintenance

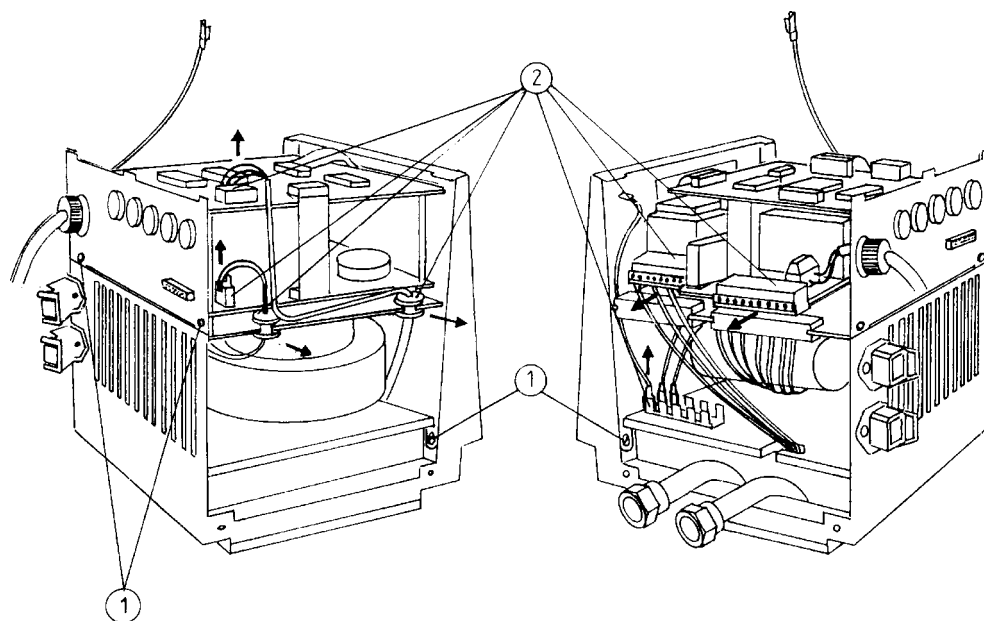
10.1 *Safety notes in case of repairs*

Always pull out the mains plug for all repair and cleaning operations! Repairs on the control unit with the cover removed must only be carried out by a qualified electrician.

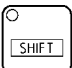
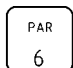
10.2 *Repair and re-initialisation*

LAUDA thermostats are largely free from maintenance. Dirty thermostatic liquid should be removed through the drain cock and replaced. If the unit should become faulty it may be advisable to return only the faulty module where appropriate.

The control unit can easily be removed after removing the cover, releasing 2 screws ① (2 turns) behind the front panel and disconnecting the electrical connections ②. The module with pump, heater, temperature probe etc. can also be separated from the bath easily.





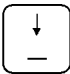

When replacing the control unit, check whether the new control unit has been programmed for the correct basic unit type. If the correct type does not appear after having switched on the mains switch, proceed as follows:

Mains switch off, press the keys  and  simultaneously and at the same time switch on the supply.

Wait until the following display appears:

L1 RK 20 K = 0 K 12 K = 1

L2 Type

Release the keys  and  and go through the menu with the key  until the required type appears. Input the code number and enter with .

The type designations are shown abbreviated, e.g. C 6 C instead of C 6 CP.

There is no provision for a fuse for the complete unit since the necessary 16A fuse corresponds to the mains fuse usually provided at the location. The control circuit of the unit has a separate fuse; a fuse 5 x 20 F4A is located in the control unit. This is accessible after removing the cover. When the fuse has blown the green lamp in the mains switch does not light up.

10.3 **Cleaning**

The unit can be cleaned using a cloth moistened with water with the addition of a few drops of (domestic) detergent. No water must find its way into the control unit. The user is responsible for any necessary decontamination if dangerous materials have been spilled on or inside the unit. This applies in particular if the unit is removed for a different use, for repair, storage etc.

The method of cleaning or decontamination is determined by the expertise of the user himself. If the user has any doubts on whether this may damage the unit he has to contact the manufacturer.

10.4 **Spares ordering**

When ordering spares please specify the equipment type and number on the label. This avoids queries and prevents the supply of wrong goods!

We shall always be happy to deal with queries, suggestions and complaints.

LAUDA DR. R. WOBSE
GMBH & CO. KG
Postfach 1251
97912 Lauda-Königshofen
Phone: (+49) (0) 9343/ 503-0
Fax: (+49) (0) 9343/ 503-222
E-mail [info @ lauda.de](mailto:info@lauda.de)
Internet <http://www.lauda.de>

11 Accessories

LAUDA Compact Temperature Thermostats C-P and K-P

LAUDA through-flow chillers

For cooling thermostats, in particular at operating temperatures below the working temperature range.

| | |
|--------|------------------|
| DLK 10 | Cat. No. LFD 105 |
| DLK 25 | Cat. No. LFD 106 |
| DLK 45 | Cat. No. LFD 107 |

Racks in stainless steel

for test tubes, centrifuge tubes etc.

Bath C 12 up to 2 racks

Bath C 20 up to 4 racks

| | | |
|---------|---|-----------------|
| RD 13 | for 56 tubes 10 - 13 dia., 80 mm immersion | Cat. No. UG 066 |
| RD 18/1 | for 33 tubes 14 - 18 dia., 80 mm immersion | Cat. No. UG 067 |
| RD 18/2 | for 33 tubes 14 - 18 dia., 110 mm immersion | Cat. No. UG 068 |
| RD 30 | for 14 tubes 24 - 30 dia., 110 mm immersion | Cat. No. UG 069 |

Bath K 20 up to 2 racks

| | | |
|---------|---|-----------------|
| RE 13 | for 56 tubes 10 - 13 dia., 80 mm immersion | Cat. No. UG 070 |
| RE 18/1 | for 33 tubes 14 - 18 dia., 80 mm immersion | Cat. No. UG 071 |
| RE 18/2 | for 33 tubes 14 - 18 dia., 110 mm immersion | Cat. No. UG 072 |
| RE 30 | for 14 tubes 24 - 30 dia., 110 mm immersion | Cat. No. UG 073 |

Bath C 6 1 rack

| | | |
|---------|---|-----------------|
| RF 18/1 | for 20 tubes 14 - 18 dia., 80 mm immersion | Cat. No. UG 074 |
| RF 18/2 | for 20 tubes 14 - 18 dia., 110 mm immersion | Cat. No. UG 075 |

Bath K 6 1 rack

| | | |
|---------|---|-----------------|
| RG 18/1 | for 20 tubes 14 - 18 dia., 80 mm immersion | Cat. No. UG 076 |
| RK 18/2 | for 20 tubes 14 - 18 dia., 110 mm immersion | Cat. No. UG 077 |

Details of other racks on request

Rising platform for retrofitting

size 250x160 mm, with continuous height adjustment suitable for:

Bath C 12 (1 platform)

Bath C 20 (2 platforms)

Bath K 20 (1 platform)

Cat. No. LCZ 012

Gable cover, stainless steel for C 20

Cat. No. LCZ 011

Level controller

for thermostating an open external bath using pressure/suction pump

Cat. No. LPZ 901

Non-return fitting

for automatic venting of the connecting hoses when thermostating an open external bath

Cat. No. UD 125

Cover plates for clear-view thermostats

| | | |
|--|----------------------|------------------|
| UD 15 K for 3 manual measuring stands | D 15 KP/ DL 15 KP | Cat. No. LTZ 017 |
| UD 20 K for 5 manual measuring stands | D 20 KP/ DL 20 KP | Cat. No. LTZ 019 |
| UD 30 K for 9 manual measuring stands | D 30 KP | Cat. No. LTZ 021 |
| UD 15 V for 2 automatic measuring stands | D 15 KP/ DL 15 KP | Cat. No. LTZ 018 |
| UD 20 V for 4 automatic measuring stands | D 20 KP/ DL 20 KP | Cat. No. LTZ 020 |
| UD 30 V for 6 automatic measuring stands | D 30 KP | Cat. No. LTZ 022 |
| UD 15 V/K for 2 manual measuring stands and 1 automatic measuring stand | D 15 KP/ DL 15 KP | Cat. No. LTZ 030 |
| UD 20 V/K for 3 manual measuring stands and 3 automatic measuring stand | D 20 KP/ DL 20 KP | Cat. No. LTZ 031 |

Cover plates made of plastic, particularly for low temperatures

| | | |
|---|----------|------------------|
| UDL 15 K for 3 manual measuring stands | DL 15 KP | Cat. No. LTZ 041 |
| UDL 15 V for 2 automatic measuring stands | DL 15 KP | Cat. No. LTZ 042 |
| UDL 20 K for 5 manual measuring stands | DL 20 KP | Cat. No. LTZ 043 |
| UDL 20 V for 4 automatic measuring stands | DL 20 KP | Cat. No. LTZ 044 |

Atherman illuminators AL

special large-area non-dazzling illuminators
free from thermal radiation

| | | |
|-------|----------------------|------------------|
| AL 15 | D 15 KP/ DL 15 KP | Cat. No. LTZ 001 |
| AL 20 | D 20 KP/ DL 20 KP | Cat. No. LTZ 002 |
| AL 30 | D 30 KP | Cat. No. LTZ 003 |

Remote operation FBC

without cable set

The electronic control unit is placed separately
from the unit and linked to it by cable.

Cable set for FBC length 5 m

Cable set for FBC length as specified

Cat. No. LCZ 960

Cat. No. UK 235

Cat. No. UK 238

Nipples

for pump connections

13 mm dia., 10 mm int. dia.

11 mm dia., 7 mm int. dia.

Cat. No. HKO 026

Cat. No. HKO 025

LAUDA Pt 100 platinum resistance thermometers
to DIN IEC 751 Class A for external control
and other temperature measurement

Pt 100-42

Cat. No. ETP 049

all-glass version with NS 14/23 ground taper DIN 12242
Temp. range -100...300°C
50% response time 0.8 sec
Overall length approx. 115 mm
4-wire circuit
Fig. 1

Pt 100-44

Cat. No. ETP 007

all-glass version with NS 14/23 ground taper DIN 12242
Temp. range -100...300°C
50% response time 0.8 sec
Overall length approx. 320mm
Fig. 2

Pt 100-66

Cat. No. ETP 008

as Pt 100-44
Overall length approx. 430mm
Fig. 2

Pt 100-90

Cat. No. ETP 050

stainless steel protection tube 4 mm dia.
Temp. range -100...300°C
50% response time 1.5 sec
Overall length approx. 120 mm
4-wire circuit
Fig. 3

Pt 100-70

Cat. No. ETP 009

stainless steel protection tube 4 mm dia.
Temp. range -200...300°C
50% response time 1.5 sec
Overall length approx. 290 mm
4-wire circuit
Fig. 3

Pt 100-92

stainless steel protection tube 4 mm dia.
with attached Silicone cable 2 m long and plug
Temp. range -100...200°C
50% response time 3 sec
Overall length approx. 250 mm
4-wire circuit
Fig. 4

Cat. No. ETP 051

Connecting cable

with 4-pin plug for external control on
all C-and K-units and for digital thermometer
for Pt 100-44 and Pt 100-66

Length 1.5m

Length as specified

for Pt 100-42. Pt 100-70. Pt 100-90

Length 1.5m

Length as specified

Cat. No. UK 048

Cat. No. UK 213

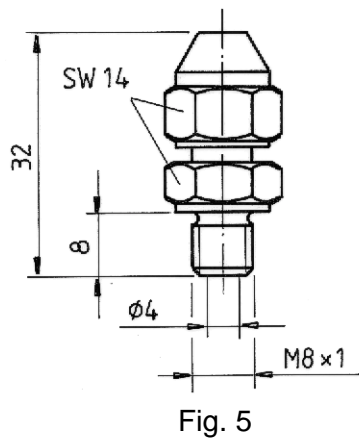
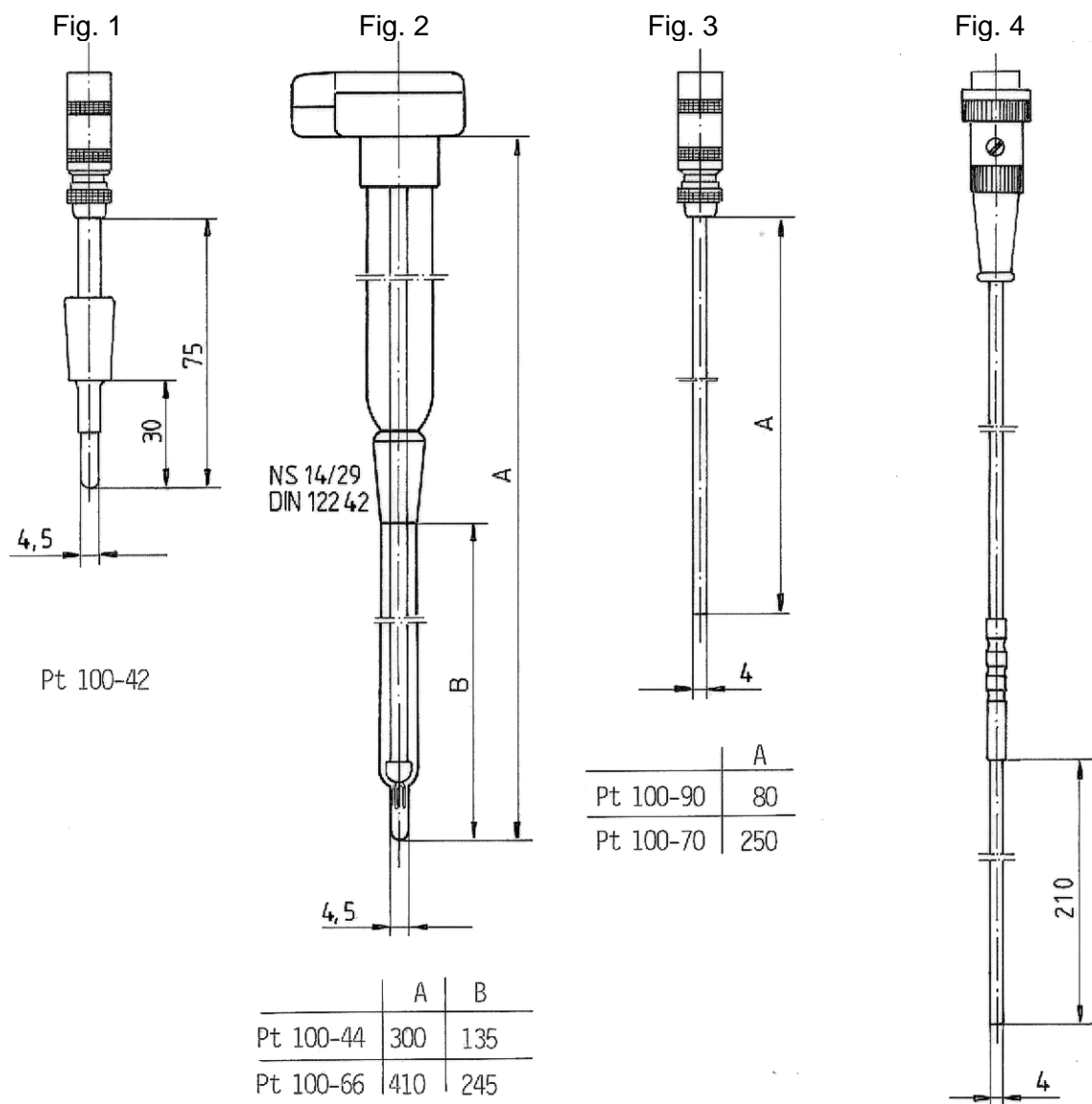
Cat. No. UK 047

Cat. No. UK 212

Screw clamp fitting

stainless steel, with Teflon pressure ring
for Pt 100 resistance thermometer 4 mm dia.
Fig. 5

Cat. No. HX 078



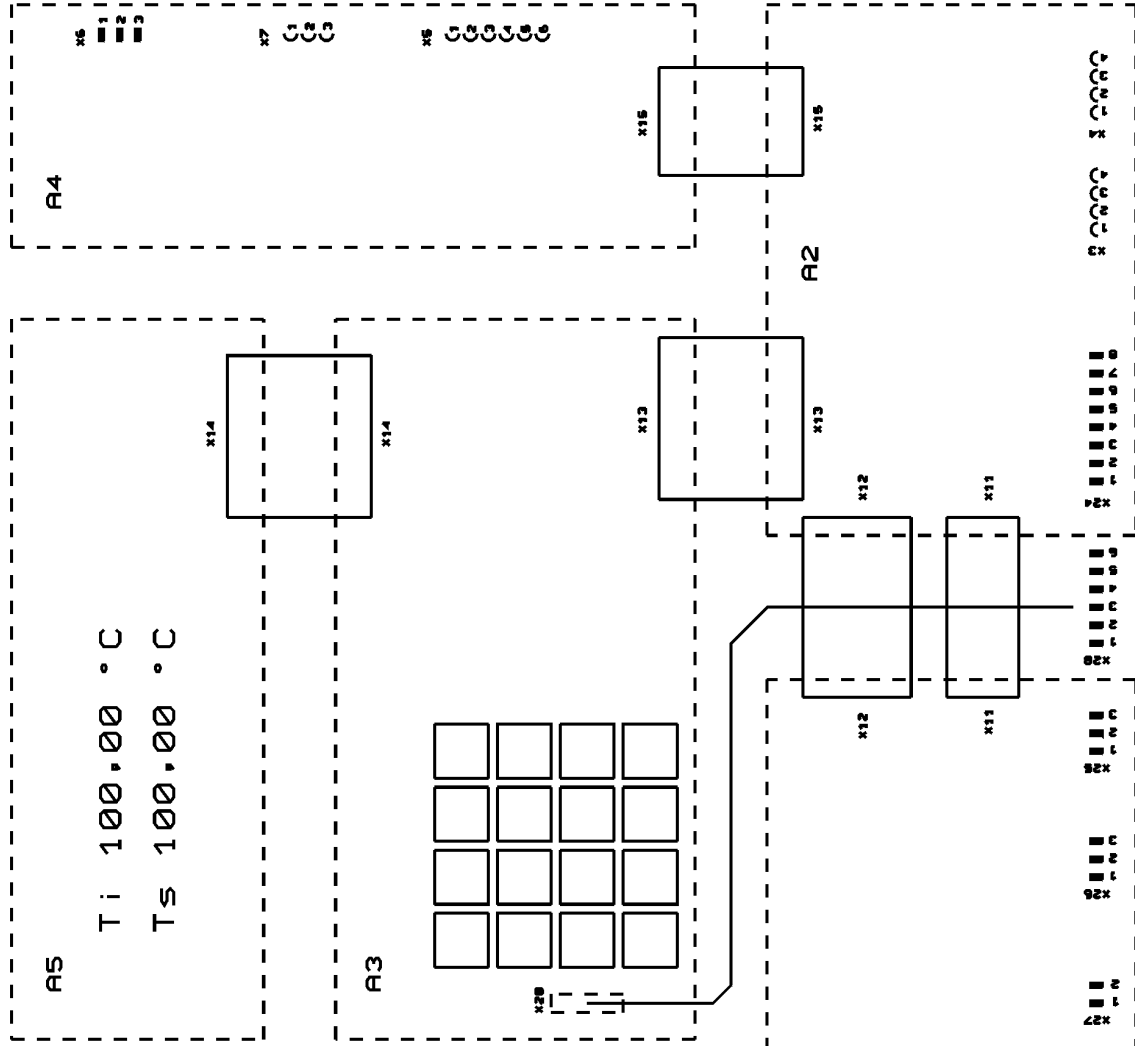
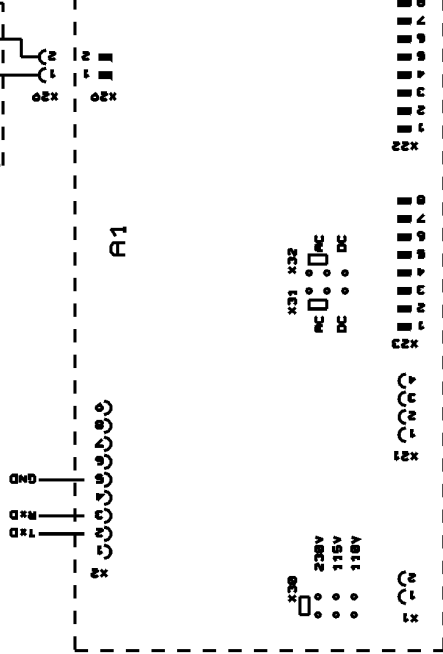
Für diese Zeichnung behalten wir uns alle nicht ausdrücklich eingeräumten Nutzungsrechte einschließlich der Vervielfältigung und Weitergabe an Dritte vor. Technische Änderungen vorbehalten

| Netzspannung | 100V | 230V |
|--------------|------|------|
| 100V | X | |
| 115V | X | |
| 230V | | X |

| x31 . x32 | AC | DC |
|-----------------|----|----|
| CxP . KxP | X | |
| DxP . DLxP | X | |
| RCxP . RKxP | X | |
| RLxP . RUL80 | | X |
| RUKxx | | X |
| R400P | X | |
| R403PL . R406PL | X | |

Anschluß

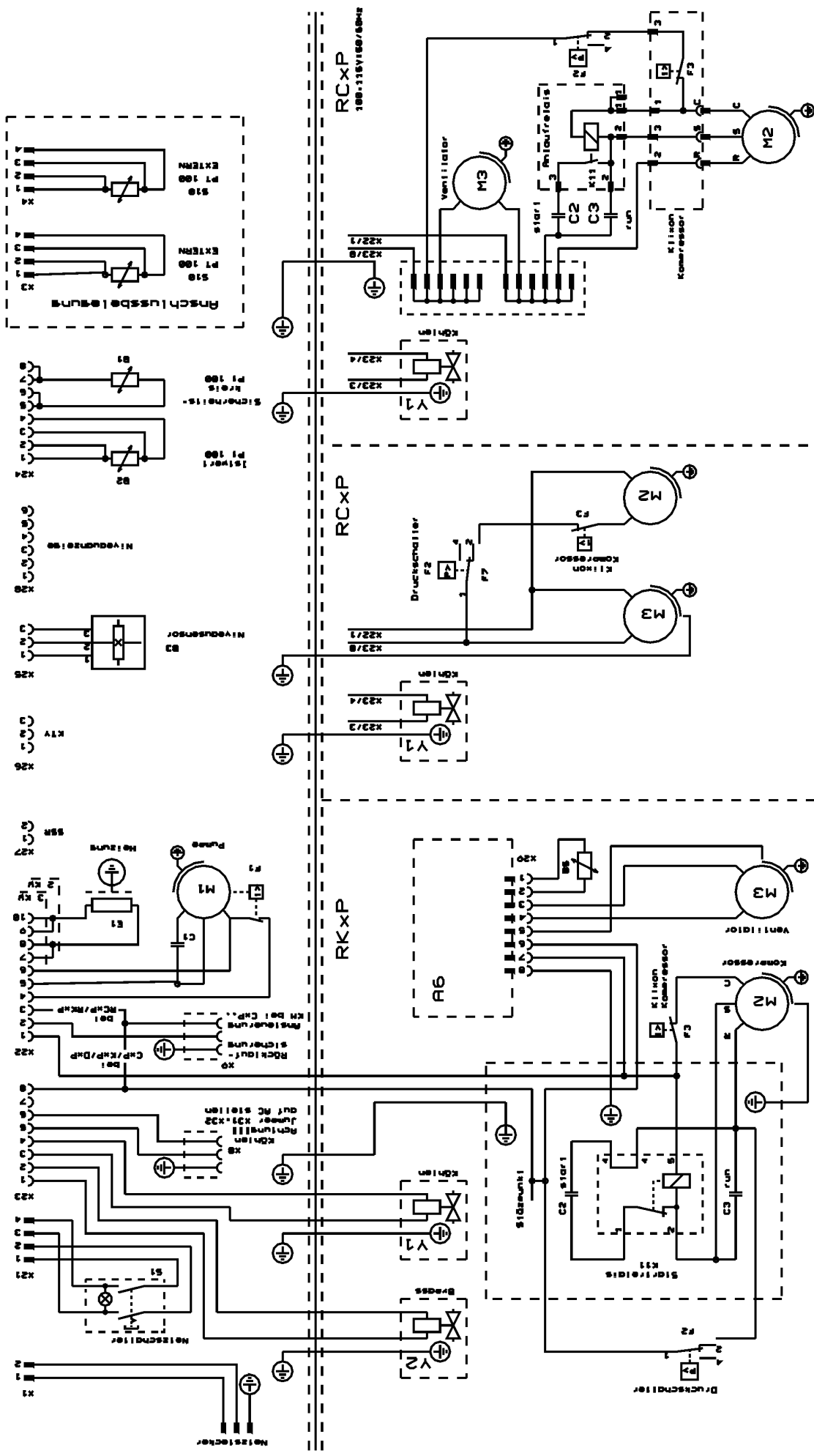
RS 2353



ab Serien. Y10...

[illegible]

Für diese Zeichnung behalten wir uns alle nicht ausdrücklich eingeräumten Nutzungsrechte einschließlich der Vervielfältigung und Weitergabe an Dritte vor. Technische Änderungen vorbehalten



| | | | |
|---|----------|----------------------------|----------|
| 14.02.1998 | | 19.03.1 | |
| Datum | Bearb. | Urspr. | Erstg. |
| 08/187 | 18.12.98 | Datum | Erstg. |
| Zust. | Benutzer | Name | Benutzer |
| ab Serie x01001 (UY122 ab Version 2.14) | | | |
| Schaltplan | | CXP/KXP/DXP/DLXP/RCXP/RKXP | |
| Schéma de connexions | | Blatt 2 | |
| Circuit diagram | | 8 | |
| LAUDA | | Dr. R. Moberg GmbH & Co KG | |

Geräteliste Schaltplan
List of parts Circuit diagram
Liste de schéma connexions
230V; 50Hz

C 6 CP
C 12 CP
C 20 CP

gültig ab Serie X01
at serial no.
à partir

| Teil-Nr. Part No. Pièce no. | Bezeichnung | Designation | Désignation | Bestell-Nr. Ref.-No. No.Ref. |
|--|---------------------------------------|---|---|---|
| A 1 | Leiterplatte „MP-Netz“ | Printed circuit board „MP-Mains“ | Circuit imprimé „MP-Secteur“ | UL 382-1B |
| A 2 | Leiterplatte „CPU“ | Printed circuit board „CPU“ | Circuit imprimé „CPU“ | UL 383-1B |
| A 3 | Leiterplatte „Anzeige-Tastatur“ | Printed circuit board „Indication-Keyboar | Circuit imprimé „Affichage-Clavier“ | UL 384-1 |
| A 4 | Leiterplatte „Analog-Ausgang“ | Printed circuit board „Analog-Exit“ | Circuit imprimé „Analogue-Sortie“ | UL 385-1 |
| A 5 | Leiterplatte „Anzeige-Display“ | Printed circuit board „Indication-Display“ | Circuit imprimé „Affichage-Display“ | EAO 091 |
| B 1 | Pt 100 Fühler Sicherheitskreis | Pt 100 Probe Safety circuit | Pt 100 Sonde Circuit sécurité | ETP 046 (US 060) |
| B 2 | Pt 100 Fühler Regelung | Pt 100 Probe Controller | Pt 100 Sonde Réglage | ETP 046 (US 060) |
| B 3 | Niveausensor | Level sensor | Niveau sensor | EKS 034 |
| C 1 | Motorkondensator | Motor condenser | Condensateur moteur | ECA 007 |
| E 1 | Heizkörper | Heater | Corps de chauffe | EH 148 |
| F 1 | Übertemperaturschutz (Umwälzpumpe) | Overtemperature protection (Circulating pump) | Protection de surpression (Pompe de circulation) | ----- |
| M 1 | Pumpenmotor | Pump motor | Moteur de pompe | EM 093 |
| S 1 | Netzschalter | Mains switch | Interrupteur secteur | EST 032 |
| X 1 | Netzanschluss / Netzkabel | Mains connection / Mains cable | Branchement secteur / Câble de secteur | EKN 008 |
| X 8 | Anschlussbuchse Kühlen | Connection socket Cooling | Douille de jonction réfroidissement | EQK 004 |
| X 9 | Anschlussbuchse Rücklaufsicherung | Connection socket Reflow security valve | Douille de jonction Protection de refoulement | EQD 037 + EQZ 006 |
| X 22 | Steckleiste 10pol. | Plug strip terminal | Réglette à fiches mâles | EZK 056 |
| X 23 | Steckleiste 8pol. (Magnetventil) | Plug strip terminal (Solenoid valve) | Réglette à fiches mâles (Vanne solenoide) | EZK 057 |

| | | | | K 6 KP K 20 KP | K 12 KP |
|-----------------------------------|---------------------------------------|---|---|------------------------------------|------------------------------------|
| Teil-Nr. Part No. Pièce no. | Bezeichnung | Designation | Désignation | Bestell-Nr. Ref.-No. No.Ref. | Bestell-Nr. Ref.-No. No.Ref. |
| A 1 | Leiterplatte „MP-Netz“ | Printed circuit board „MP-Mains“ | Circuit imprimé „MP-Secteur“ | UL 382-1B | UL 382-1B |
| A 2 | Leiterplatte „CPU“ | Printed circuit board „CPU“ | Circuit imprimé „CPU“ | UL 383-1B | UL 383-1B |
| A 3 | Leiterplatte „Anzeige-Tastatur“ | Printed circuit board „Indication-Keyboar | Circuit imprimé „Affichage-Clavier“ | UL 384-1 | UL 384-1 |
| A 4 | Leiterplatte „Analog-Ausgang“ | Printed circuit board „Analog-Exit“ | Circuit imprimé „Analogue-Sortie“ | UL 385-1 | UL 385-1 |
| A 5 | Leiterplatte „Anzeige-Display“ | Printed circuit board „Indication-Display“ | Circuit imprimé „Affichage-Display“ | EAO 091 | EAO 091 |
| B 1 | Pt 100 Fühler Sicherheitskreis | Pt 100 Probe Safety circuit | Pt 100 Sonde Circuit sécurité | ETP 046 (US 060) | ETP 048 (US 061) |
| B 2 | Pt 100 Fühler Regelung | Pt 100 Probe Controller | Pt 100 Sonde Réglage | ETP 046 (US 060) | ETP 048 (US 061) |
| B 3 | Niveausensor | Level sensor | Niveau sensor | EKS 034 | EKS 034 |
| C 1 | Motorkondensator | Motor condenser | Condensateur moteur | ECA 007 | ECA 007 |
| E 1 | Heizkörper | Heater | Corps de chauffe | EH 150 | EH 154 |
| F 1 | Übertemperaturschutz (Umwälzpumpe) | Overtemperature protection (Circulating pump) | Protection de surpression (Pompe de circulation) | ----- | ----- |
| M 1 | Pumpenmotor | Pump motor | Moteur de pompe | EM 094 | EM 096 |
| S 1 | Netzschalter | Mains switch | Interrupteur secteur | EST 032 | EST 032 |
| X 1 | Netzanschluss / Netzkabel | Mains connection / Mains cable | Branchement secteur / Câble de secteur | EKN 008 | EKN 008 |
| X 8 | Anschlussbuchse Kühlen | Connection socket Cooling | Douille de jonction réfroidissement | EQK 004 | EQK 004 |
| X 9 | Anschlussbuchse Rücklaufsicherung | Connection socket Reflow security valve | Douille de jonction Protection de refoulement | EQD 037 + EQZ 006 | EQD 037 + EQZ 006 |
| X 22 | Steckleiste 10pol. | Plug strip terminal | Réglette à fiches mâles | EZK 056 | EZK 056 |
| X 23 | Steckleiste 8pol. (Magnetventil) | Plug strip terminal (Solenoid valve) | Réglette à fiches mâles (Vanne solenoide) | EZK 057 | EZK 057 |

Geräteliste Schaltplan
List of parts Circuit diagram
Liste de schéma connexions
230V; 50Hz / 230V; 60Hz

D 15 KP
D 20 KP
D 30 KP

DL 15 KP
DL 20 KP

gültig ab Serie X01
at serial no.
à partir

| Teil-Nr. Part No. Pièce no. | Bezeichnung | Designation | Désignation | Bestell-Nr. Ref.-No. No.Ref. |
|--|---------------------------------------|---|---|---|
| A 1 | Leiterplatte „MP-Netz“ | Printed circuit board „MP-Mains“ | Circuit imprimé „MP-Secteur“ | UL 382-1B |
| A 2 | Leiterplatte „CPU“ | Printed circuit board „CPU“ | Circuit imprimé „CPU“ | UL 383-1B |
| A 3 | Leiterplatte „Anzeige-Tastatur“ | Printed circuit board „Indication-Keyboar | Circuit imprimé „Affichage-Clavier“ | UL 384-1 |
| A 4 | Leiterplatte „Analog-Ausgang“ | Printed circuit board „Analog-Exit“ | Circuit imprimé „Analogue-Sortie“ | UL 385-1 |
| A 5 | Leiterplatte „Anzeige-Display“ | Printed circuit board „Indication-Display“ | Circuit imprimé „Affichage-Display“ | EAO 091 |
| B 1 | Pt 100 Fühler Sicherheitskreis | Pt 100 Probe Safety circuit | Pt 100 Sonde Circuit sécurité | ETP 048 (US 061) |
| B 2 | Pt 100 Fühler Regelung | Pt 100 Probe Controller | Pt 100 Sonde Réglage | ETP 048 (US 061) |
| B 3 | Niveausensor | Level sensor | Niveau sensor | EKS 034 |
| C 1 | Motorkondensator | Motor condenser | Condensateur moteur | ECA 007 |
| E 1 | Heizkörper | Heater | Corps de chauffe | EH 154 |
| F 1 | Übertemperaturschutz (Umwälzpumpe) | Overtemperature protection (Circulating pump) | Protection de surpression (Pompe de circulation) | ----- |
| M 1 | Pumpenmotor | Pump motor | Moteur de pompe | EM 096 |
| S 1 | Netzschalter | Mains switch | Interrupteur secteur | EST 032 |
| X 1 | Netzanschluss / Netzkabel | Mains connection / Mains cable | Branchement secteur / Câble de secteur | EKN 008 |
| X 8 | Anschlussbuchse Kühlen | Connection socket Cooling | Douille de jonction réfroidissement | EQK 004 |
| X 9 | Anschlussbuchse Rücklaufsicherung | Connection socket Reflow security valve | Douille de jonction Protection de refoulement | EQD 037 + EQZ 006 |

Each repair requiring the opening of the control part by means of tools and each work at the electronic part may only be done by a trained technician !

| Malfunction | Fault | Reason | Remedy |
|--|---|--|---|
| Green signal lamp of mains switch not alight | Control fuse has actuated | | Replace fuse on PCB "MP Mains" 5 x 20; F4A |
| | | Overload on PCB | Replace PCB "MP Mains" |
| Display: "Level too low " | Bath level too low | Evaporation; external consuming device not refilled | Fill in bath liquid; press reset button twice |
| | | Leakage in hose connections | Check hoses and their connection; if necessary replace them; fill in bath liquid; press reset button twice |
| | Level detector | | Check level detector, plugs and hall sensors; if necessary replace them; check their function carefully For USH 400(/6) fill in bath liquid up to the level mark |
| Display: " Temperature too high " | Temperature probe | Temperature difference between two probes > 15°C | Replace bath temperature probe (double Pt 100). For USH 400(/6) please check both single Pt 100. |
| | Temperature of bath liquid above over-temperature switch-off point (To) | Overtemperature switch-off point (To) set too low | Press reset button; adjust overtemp. switch-off point (To) at an higher value; press reset button |
| | | Triac or triggering | Replace triac or PCB "MP Mains" |
| Display: " External fault " | Contact input 14 N "Fault " is used | Pins ½ not connected | Reason for malfunction caused by external system |
| | | No signal transmitter connected to socket 14 N, though selection at PAR level is activated | Switch off function " Contact input Fault " at PAR level; see operating instructions |
| Display: " Ts > To | Wrong inputs | Setpoint adjusted above overtemperature switch-off point (To) | Adjust overtemperature switch-off point (To) at an higher temperature; pay attention to bath liquid, flash point etc. !! |
| | | Overtemperature switch-off point adjusted below setpoint (Ts) | First adjust setpoint (Ts) at a lower value then set requested overtemperature switch-off point |
| Display: "Ts < Tu" | Wrong input | Setpoint adjusted below low-temperature switch-off point (Tu) | Adjust low-temperature switch-off point (Tu) at a lower value |
| | | Low-temperature switch-off point adjusted above setpoint (Ts) | First adjust setpoint at an higher temperature then set requested low-temperature switch-off point (Tu) |

| Malfunction | Fault | Reason | Remedy |
|--|--------------------|---|---|
| Setpoint Ts is not adjustable but will disappear | Operating error | Setpoint is determined by the analogue input; see right side of display L2: A | Switch off analogue input |
| Sound signal appears when a button is pressed | | Another function blocks the keyboard e.g.: programme runs; RS 232 active; Parameter etc. | Leave the function or press ESC (RS 232) (R appears on the right side of display L2); Stop access and press the button "Local" |
| Display: "Tu – Cursor flashes" acoustic signal switches compressor off after 1 min. | Wrong input | Actual value is \leq Tu, resp. setpoint adjusted too close to low-temperature switch-off point Tu; bath temperature (Ti) falls below Tu | Set low-temperature switch-off point (Tu) at a lower value |
| Display "Out of range" | Wrong input | Tried to enter values being out of admissible ranges; Ts, To, Tu being out of operating temp. | Choose the right values taking into consideration their limitations; check bath liquid or configuration after having switched on the unit |
| | | Programme input out of operating temperature range of the unit | Enter admissible values |
| | | Value for Xp, Tn, Tv above 199,9 | Enter admissible values |
| Tu or To is not adjustable; "Out of range " | | Input values are outside of temperature limits of the type of unit or initialization does not fit to the type of unit | Reinitialize type of unit: (see operating instructions "Maintenance"); may also return Default-values |
| Display: "TA" (only for USH 400(/6) | | Motor chamber temperature > 55 °C | Surrounding temperature of the part of the thermostat may be too high (see operating instructions 5.2.8) |
| Display: ↓ (only for USH 400(/6) RUL and RUK) | | Level too low | Fill in bath liquid ; (see operating instructions 5.2.3) |
| Display: ↑ (only for USH 400(/6) RUL and RUK) | | Level in the vessel is close to overflow; heater switches off | Either reduce the amount of thermal liquid or install an other vessel Attention: HOT !!! |
| Display: "Upper limitation of oil > limitation of unit" (cancelled beginning with software version 1.06) | | Admissible working temperature range of bath liquid exceeds operating temperature range of the unit | Bath liquid is accepted, no other steps necessary; unit limitations valid |
| Setpoint cannot be selected by means of keyboard | | Check setpoint selection; see right side of display L2; P=Programme; A=Analogue; R=RS 232 | Switch setpoint selection to I=Internal |
| Display: "External fault – clock stop!" | Clock does not run | RAM defective | Unit has to be switched on once again; set date and time once again; see operating instructions if necessary replace RAM |

| Malfunction | Fault | Reason | Remedy |
|--|---|---|---|
| Display: "Internal Pt 100 defective" | Double Pt 100 for bath temperature or safety temperature | Interruption, short circuit or temperature deviation of bath temperature probes too important | Replace double Pt 100 for bath temperature For USH 400(/6) please check both single Pt 100 |
| Display: "Ext. Pt 100 not connected" | | Tried to switch over to external control without connecting an external Pt 100 Tried to calibrate Pt 100 being not connected | Continue to work with internal control or connect external Pt 100 for T1 or T2; look at display for control variables I, 1 or 2; check display for T1 or T2 |
| Display: "Correcting value too high" | Important deviation of Pt 100 from standard values | Value input differing from the basic value that is indicated by more than 5 K | Check temperature reference thermometer, check Pt 100; replace PCB "CPU" |
| Scale of analogue outputs shows discrepancies | | Outputs not correctly calibrated | Calibrate the analogue output channels (see operating instructions) |
| Unit does not heat though heating is indicated | Triac | | Replace triac |
| | Heater | Defective | Replace heater |
| | | Interruption | Eliminate |
| | | Electronics | Replace PCB "MP Mains" |
| | | Controller output limitation at PAR level too small | Enter higher values at PAR level (e.g. 100 %) |
| Pump does not run | Temperature safety cut-out in pump has actuated Pump stops | Motor blocked | Turn propeller of motor; if necessary replace it; clean pump |
| | | Viscosity of bath liquid too high | Use other bath liquid; wait until motor has cooled down |
| Bath temperature rises clearly above adjusted setpoint (Ts) Heating indication ON | | Controller | Replace PCB "MP Mains" or PCB "CPU " |
| Heating indication OFF | | Triac | Replace triac |
| Temperature rises slowly above adjusted setpoint Heating indication OFF | Cooling not sufficient | Heat emission of pump | Connect cooling water supply or other kind of cooling |
| Display shows wrong temperature (Ti, T1, T2) | | Temperature probe | Replace double Pt 100 or external Pt 100 T1, T2 |
| Display is dark | Temperature cut-out in transformer has actuated | Overload caused by short – circuit | Replace PCBs or control unit |

| Malfunction | Fault | Reason | Remedy |
|---|---------------------------------------|--|---|
| Unit does not work at adjusted setpoint | | Wrong control parameters adjusted | Enter new values or start autoadaptation |
| Temperature control by means of external controller not stable | | Thermal contact of bath liquid and external measuring point not sufficient | Improve circulation through external consuming device or thermal contact to ext. Pt 100 |
| Unit does not cool down | | Controller output limitation at PAR level too small | Adjust higher value at PAR level (e.g. -100 %) |
| | Compressor defective | | Replace refrigeration unit – by refrigeration engineer !! |
| | Leakage in refrigeration system | | Clear leakage, fill in refrigerant – by refrigeration engineer !! |
| | Compressor does not run | Compressor without tension | Eliminate line interruption |
| | | Triggering defective | Replace PCB "MP Mains " |
| Unit does not cool down | Compressor does not run | Pressure switch has actuated | Unscrew grid, clean condenser, blow through compressed air, improve ventilation |
| | Solenoid valves do not work correctly | Triggering defective | Replace PCB " MP Mains " |
| Compressor switches ON and OFF in regular periods; temperature constancy very bad | Condensation pressure too high | Fan defective | Replace fan motor |
| | | Fan speed (only RK) | Check speed or speed controller |
| | Condenser dirty | Dust | Unscrew grid, from the back side blow compressed air or nitrogen through condenser |
| | | Ventilation disturbed | Enlarge distance to nearby units or walls |
| | | Ambient temperature too high | Air the room |
| Insufficient cooling in the lower temperature range | | Bath liquid contaminated by condensate | Replace bath liquid to suit bath temperature |
| | Cools down to approx. 0°C only | Bath liquid not suitable (water) | Use water/ glycol mixture |

An / To / A:

LAUDA Dr. R. Wobser • LAUDA Service Center • Fax: +49 (0) 9343 - 503-222

Von / From / De :

Firma / Company / Entreprise: _____

Straße / Street / Rue: _____

Ort / City / Ville: _____

Tel.: _____

Fax: _____

Betreiber / Responsible person / Personne responsable: _____

Hiermit bestätigen wir, daß nachfolgend aufgeführtes LAUDA-Gerät (Daten vom Typenschild):

We herewith confirm that the following LAUDA-equipment (see label):

Par la présente nous confirmons que l'appareil LAUDA (voir plaque signalétique):

| Typ / Type / Type : | Serien-Nr. / Serial no. / No. de série: |
|---------------------|---|
| | |

mit folgendem Medium betrieben wurde

was used with the below mentioned media

a été utilisé avec le liquide suivant

Darüber hinaus bestätigen wir, daß das oben aufgeführte Gerät sorgfältig gereinigt wurde, die Anschlüsse verschlossen sind, und sich weder giftige, aggressive, radioaktive noch andere gefährliche Medien in dem Gerät befinden.

Additionally we confirm that the above mentioned equipment has been cleaned, that all connectors are closed and that there are no poisonous, aggressive, radioactive or other dangerous media inside the equipment.

D'autre part, nous confirmons que l'appareil mentionné ci-dessus a été nettoyé correctement, que les tubulures sont fermées et qu'il n'y a aucun produit toxique, agressif, radioactif ou autre produit nocif ou dangereux dans la cuve.

| Stempel Seal / Cachet. | Datum Date / Date | Betreiber Responsible person / Personne responsable |
|---------------------------|----------------------|--|
| | | |